

KERN & SOHN GmbH
Balances, Test weights, Microscopes,
DAkkS calibration laboratory
Ziegelei 1
72336 Balingen
Germany

Tel. +49 7433 9933-0
Fax +49 7433 9933-149

info@kern-sohn.com
www.kern-sohn.com

Printed in Germany
by KERN & SOHN GmbH
zcm-en-kp-20211



MEDICAL SCALES

for clinics, doctors' surgeries, care facilities and vets



1769
Ancestor Johann Jakob Sauter built the Hahn inclination scale out of iron, a foundation stone for weighing for the balance industry in Southern Germany



1844
KERN is founded – precision balances are produced



1863
A proud Gottlieb Kern with his staff



1880
Pharmaceutical balance with Aesculap



1923
Inflation – KERN wages are paid with self printed currency



1980
The electronic balance ousts mechanical devices



1994
Accredited DKD laboratory (ISO 17025)



2000
New premises in Balingen



2002
Existing QM system certification in accordance with DIN EN ISO 9001:2000 standards



2007
Approval for the manufacture of medical products (DIN EN 13485 and 93/42/EEC)



2008
Authorisation for initial verification by the manufacturer (2009/23/EC)



2009
Approval for the manufacture and sale of height rods (DIN EN 13485 and 93/42/EEC)



2012
Verification point for non-automatic balances and test weights.

New customer portal
www.kern-sohn.com goes live



2014
Expansion of the product range to include optical instruments (microscopes and refractometers)



2015
Inauguration of Ziegelei 2.0 with computer-controlled high-bay warehouse



2017
Come with KERN into the digital future: Expansion of the model ranges compatible with Industry 4.0, as well as the related services



2019
Significant anniversary year! 25 years of accredited DAkkS/DKD laboratory 175 years of KERN & SOHN 250 years of balance manufacture in the Sauter family-owned company



2020
Construction of Ziegelei 3.0, extension of administration building

**NEW WEBSHOP
ONLINE!**
www.kern-sohn.com

PROFESSIONAL MEASURING

2021



Legal provisions concerning medical products

Legal notice

Medical personal scales are understood to be all non-automatic scales that are intended to be used to determine the body weight in medical applications (e.g. patient, baby, chair and bed scales, verified according to class). These scales are regulated by European Directive 2014/31/EU.

- Medical personal scales are also governed by the regulations for medical devices that are placed on the market in accordance with European Directive 93/42/EEC (for Germany: Medical Products Law [MPG])

Types of scales for medical applications that require mandatory verification:

APPROVED MEDICAL DEVICE

- Personal scales (incl. chair and wheelchair scales) in hospitals
- Baby scales and mechanical birth weight scales
- Bed scales
- Scales in dialysis units

Verified medical personal scales with a maximum resolution of 3000 digits are classified as Verification Class III (see type plate).

The re-verification deadlines depend on the location of use and are given in the Verification Ordinance. The operating manual specifies the re-verification deadlines for Germany at the time of going to press. The operator of the scale is obliged to ensure that the verification stamps are not damaged and that the statutory re-verification deadlines are observed.

KERN & SOHN GmbH does not accept responsibility for this.

Scales not specifically intended as a medical device and without verification:

These scales do not have an approval according to Directive 93/42/EEC (or the German Medical Products Law). They are thus not permitted for use in medical applications. They are used for self-checks by the patient and are not verified.

A mandatory verification is not required for the following intended uses:

PROFESSIONAL CARE

- Body weight scales in pathology (not a medical application according to § 7b Paragraph 2 No. 4 of the Verification Ordinance)
- Personal and baby scales (the latter can usually be borrowed in pharmacies)
- Baby scales for midwives
- Personal scales during blood collection for the preparation of blood products

(Information from: "Die Eichverwaltung informiert, Waagen in der Heilkunde" (Information from the Verification Office: Scales in medical applications) dated 19.04.2013)

HOME CARE

Non-verified scales are not permitted to be used in medical applications. They are used for self-monitoring of the health status. Examples of these include personal scales, body fat scales, etc. Such scales can be used to determine weight trends so that a doctor can be consulted in case of doubt. The patient is then weighed with verified scales at the medical practice to enable diagnosis.

Other information

Balances with an EC Type-Approval which have been distributed without verification, cannot be subsequently verified. Non-verified balances are not acknowledged as a medical product in accordance to norm 93/42/EEC and therefore, may not be used in medical applications.

Gravitational force

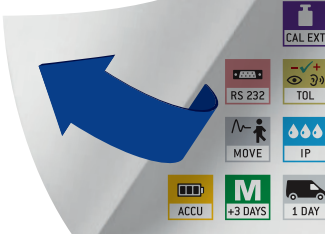
Electronic balances are subject to location regulations due to the earth's gravitational force. If the location of use is changed, the electronic balance must be re-verified and re-calibrated.

Übersicht über die Symbole
El GEÄIÄHÄ€Mä@äO{ ä&ä&ä








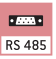













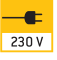











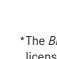

For an overview of the symbols,
see the reverse side

Images, page:
6 ©Syda Productions/shutterstock.com, 8 ©Africa Studio/shutterstock.com, 12 ©Suteren/shutterstock.com,
14 ©megaflopp/shutterstock.com, 16 ©Goran Bogicevic/shutterstock.com, 18 ©Africa Studio/shutterstock.com,
20 ©Idutko/shutterstock.com, 22 ©didesign021/shutterstock.com, 28 ©Rido/shutterstock.com,
32 ©Syda Productions/shutterstock.com, 34 ©LightField Studios/shutterstock.com, 36 ©Robert Kneschke/shutterstock.com,
38 @ JP WALLET/shutterstock.com, 40 ©New Africa/shutterstock.com, 42 ©ALPA PROD/shutterstock.com,
47 @ YAKOBCHUK VIACHESLAV/shutterstock.com, 53 ©Africa Studio/shutterstock.com

Übersicht über die Symbole
El GEÄIÄHÄ€Mä@äO{ ä&ä&ä



KERN Pictograms

 Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required	 ZERO: Resets the display to "0"	 Monocular Microscope: For the inspection with one eye
 Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	 Hold function: When patients do not stand, sit or lie completely still, a stable weight is calculated using an average weight	 Binocular Microscope: For the inspection with both eyes
 Data interface RS-232: To connect the balance to a printer, PC or network	 Hold function: When the weighing conditions are unstable, a stable weight is calculated as an average value	 Trinocular Microscope: For the inspection with both eyes and the additional option for the connection of a camera
 RS-485 data interface: To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible	 Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram	 Abbe Condenser: With high numerical aperture for the concentration and the focusing of light
 Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals	 Suspended weighing: Load support with hook on the underside of the balance	 Halogen illumination: For pictures bright and rich in contrast
 WIFI data interface: To transfer data from the balance to a printer, PC or other peripherals	 Battery operation: Ready for battery operation. The battery type is specified for each device	 LED illumination: Cold, energy-saving and especially long-life illumination
 Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.	 Rechargeable battery pack: Rechargeable set	 Fluorescence illumination for compound microscopes: With 100 W mercury lamp and filter
 Statistics: using the saved values, the device calculates statistical data, such as average value, standard deviation etc.	 Battery operation rechargeable Prepared for a rechargeable battery operation	 Fluorescence illumination for compound microscopes: With 3 W LED illumination and filter
 PC Software: to transfer the measurements from the device to a PC	 Universal mains adapter: with universal input and optional input socket adapters for A) EU, CH B) EU, CH, GB, USA	 Phase contrast unit: For a higher contrast
 GLP/ISO-Protokoll: With date and time. Only with KERN printers	 Mains adapter: 230V/50Hz in standard version for EU. On request GB, AUS or USA version available	 Darkfield condenser/unit: For a higher contrast due to indirect illumination
 KERN Communication Protocol (KCP): It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems	 Power supply: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request	 Polarising unit: To polarise the light
 Piece counting: Reference quantities selectable. Display can be switched from piece to weight	 Weighing principle: Strain gauges Electrical resistor on an elastic deforming body	 Infinity system: Infinity corrected optical system
 Totalising level A: The weights of similar items can be added together and the total can be printed out	 Peak hold function: capturing a peak value within a measuring process	 Automatic temperature compensation: For measurements between 10 °C and 30 °C
 Weighing units: Can be switched to e.g. nonmetric units at the touch of a key. Please refer to website for more details	 Push and Pull: the measuring device can capture tension and compression forces	 Verification possible: The time required for verification is specified in the pictogram
 Weighing with tolerance range (Check weighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model	 Integrated scale: In the eyepiece	 Package shipment: The time required for internal shipping preparations is shown in days in the pictogram
	 360° rotatable microscope head	 Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.

KERN – Your partner for Medical scales

Quality & Experience

We are absolutely committed to making products available to you which are of excellent built quality, are extremely robust, have a long service life and are easy to use too. These three factors: quality, durability and ease of operation have been hallmarks of our products for over 175 years.

Conformity

Particularly in the medical environment, balances are subject to extensive standards and regulations, which serve to protect the users and patients. Over and above these standards, our certified quality management system, as well as our 20 trained medical product advisers ensure rapid and competent product advice on all aspects of medical weighing technology.

Your advantages

fast

- 24 hour dispatch service for products in stock – order today, on its way tomorrow
- Sales & service hotline from 8:00 am to 6:00 pm


reliable


- Up to 3 years warranty
- Precision in weighing technology for more than 175 years

 **Order hotline**
+49 7433 9933-0

 **Service hotline**
+49 7433 9933-199

 **Calibration hotline**
+49 7433 9933-196

 **Our team of consultants will assist you**
from Monday to Friday
from 8:00 am to 6:00 pm

 **Online shop**
www.kern-sohn.com

 **E-mail order**
info@kern-sohn.com

 **Fax order**
+49 7433 9933-146

 **Our team of consultants will assist you**
from Monday to Friday
from 8:00 am to 6:00 pm

Promptness

The right product, at the right time, in the right place – that is our strength. With our sophisticated internal logistics system we can make sure your order is ready for delivery within 24 hours (plus calibration time if applicable). So you can react quickly to your customers' requirements and requests.

Professionalism and trust

On the subject of competent advice: Competence coupled with reliability – for centuries this has been the guiding principle of our company. This is why KERN sells its products exclusively through authorised dealers.

Albert Sauter, Managing Director
[Balingen, December 2020]

Retailer information

Sales conditions

All prices are valid as of January 1st 2021 until a new version of the KERN catalogue is released. In Europe, all prices do not include the applicable V.A.T.

Extract from general terms and conditions: Court of jurisdiction/ Legal domicile: 72336 Balingen, Germany; Commercial register N°: HRB 400865, AG Stuttgart; Managing director: Albert Sauter. For the full Terms and Conditions, please refer to the website. www.kern-sohn.com/en/kern/agbs.html

Price changes and product changes are likely in individual cases, due to product modifications as well as error.

Sale or return: within 14 days of purchase. Not valid for customised items, such as special orders or test services, such as e.g. calibration, verification etc. Depending on the time and effort involved, there may be processing and storage costs, please ask for details.

Warranty:
3 years for products with list price ≥ € 400.00
2 years for components as well as products with list price < € 400.00
Does not apply to consumables such as batteries, rechargeable battery packs, etc.

Balances in the legally regulated range (verification) must comply with the commercial error limits during the verification validity period – these usually amount to twice the permissible error limits. If this verification period expires, a reverification must be carried out. If an adjustment of the scale to comply with the permissible error limits is necessary to pass this reverification, this does not constitute a warranty claim.

Delivery Conditions: we supply ex works Balingen, e.g. the transport costs are invoiced. Delivery is usually via courier service, or by truck. Please ask for prices. Any goods supplied, remain KERN's property until complete payment for the goods sold has been received.

All dimensions are approximate. Subject to technical changes without prior notice as a result.

Measuring technology, microscopy and testing services from a single source

				
Medical scales catalogue	Balances & Test service catalogue	Microscopes & refractometers catalogue	SAUTER measuring equipment catalogue	KERN DAKkS calibration service brochure

KERN Medical Scales 2021









Height measuring systems

MS... · MBA	Tape for measuring circumference, Height rods	04, 05
-------------	---	--------

Baby scales

MBC		Baby scale	06
MBA		Baby scale	08
MBD		Baby scale	10

Personal floor scales

MPE · MPE-E	 	Personal floor scale _____	12
MPS	 	Personal floor scale _____	14
MPB		Personal floor scale _____	16
MPC	 	Personal floor scale/Obesity scale _____	18
MPD · MPD-E	 	Personal floor scales/Doctor scales _____	20
MPK · MPL		Personal floor scales _____	22
MPI		Personal floor scale _____	24
MGC		Mechanical personal floor scale _____	25




Bathroom scales

MGD	Bathroom scale	26
-----	----------------	----

Body analysis scales

MFB Body analysis scale _____ 27

Wheelchair platform scales

MWA		Wheelchair platform scale	28
MWB		Wheelchair platform scale	31
MWS		Wheelchair platform scale	32

Stretcher scales

MWS-L		Stretcher scale	32
-------	---	-----------------	----

Handrail scales

MTA Handrail scale _____ 34

Chair scales

MCB		Chair scale	36
MCC		Chair scale	38
MCD		Chair scale	40

Laundry cart scales

NFN	Laundry cart scale	43
-----	--------------------	----

Organ scales

FXN Organ scale _____ 44

Dynamometers

MAP	Digital dynamometer_____	45
FK · FH	Digital dynamometers for medical applications_____	46

Microscopes

OBE-12 · 13	Compound microscope	48
OBF-1 · OBL-1	Compound microscopes	49
OBN-1	Compound microscope	50
OCM-1	Inverted microscopes	51

Refractometers

ORA	Analogue refractometer	52
-----	------------------------	----

Veterinary scales

HCB · HCN	Hanging scales_____	54
EMB	Small pet scale_____	55
FOB	Stainless steel small pet scales_____	56
EOB · EOE	Veterinary scales _____	57
EOS	Veterinary scale _____	58

KERN Models A-Z

E

EMB_____	55
EOB/EOE_____	57
EOS_____	58

F

FH/FK _____	46
FOB _____	56
FXN _____	44

H

HCB/HCN	54
---------	----

M

MAP	45
MBA-A01	04
MBC	06
MBD	10
MCB	36
MCC	38
MCD	40
MGC	25
MGD	26
MFB	27
MPB	16
MPC	18
MPD/MPD-E	20
MPE/MPE-E	12
MPI	24
MPK/MPL	22
MPS	14
MSB	05
MSC	05
MSF	05
MSW	04
MTA	34
MWA	28
MWB	31
MWS	32
MWS-L	32

N


NFN	43
-----	----

0

OBE	48
OBF/OBL	49
OBN	50
OCM	51
ORA	52

■ Approved Medical Device

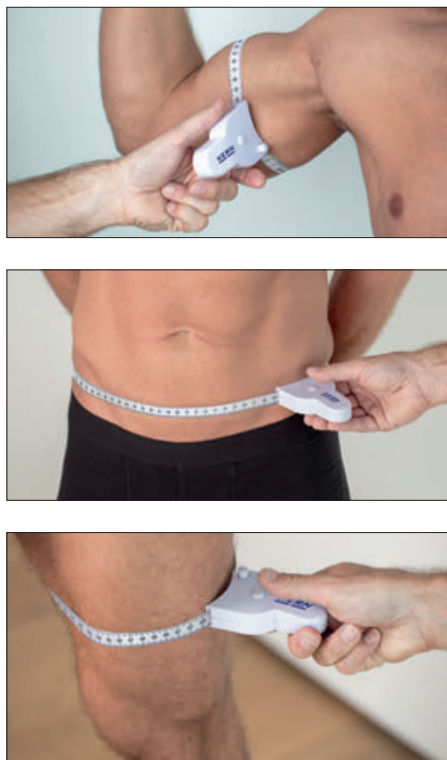- EC type approval



Circumference measuring tape and height rods for professional use in medical diagnostics

1 KERN MSW 200

Convenient tape for measuring circumference with convenient pull-out mechanism

Features

- Continuous pull-out measuring tape for e.g. precise measurement around the body
- Ideal to determine the WHR (Waist-to-Hip-Ratio), which can give information about the distribution of the body's fat deposits using an indirect determination of abdominal fat
- Useful diagnosis instrument
- With integrated feed mechanism
- The easiest handling:
 - Pull out the measuring tape,
 - Hang it in the recess provided for this,
 - Operate the feed mechanism,
 - Read measurement value with millimetre precision
- Measuring tape length 2,05 m (80 inch)
- Net weight approx. 0,05 kg
- **Note:** The models are only delivered in a set of 5 units. The price refers to a delivery of 5 items. Cannot be delivered individually.
- KERN MSW 200S05, € 31,-

2 KERN MBA-A01

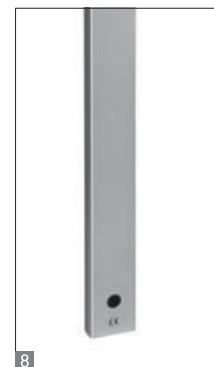
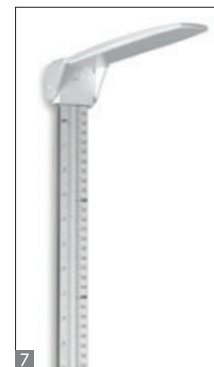
Portable mechanical height rod

Features

- Robust construction
- Compact size
- Easy and hygienic cleaning
- For babies up to a maximum of 80 cm
- Measuring range 30–80 cm (11,8–31,5 inch); increments of 1 mm (0,04 inch)
- Readability on scale with moveable stop
- Large guide surfaces make handling easier (aligning, shifting, reading)
- Overall dimensions max. W×D×H 850×240×60 mm
- Net weight approx. 0,36 kg
- KERN MBA-A01, € 110,-

STANDARD





When every centimetre counts – mechanical height rods approved as a medical device for professional use in medical diagnostics

3 KERN MSB 80

Portable mechanical height rod

Features

- Approved as a medical device according to 93/42/EEC
- For babies of up to 80 cm
- Measuring range 10–80 cm; increments of 1 mm
- Large guide surfaces make handling easier (aligning, shifting, reading)
- 4 Measuring scale begins with 0 at right and left, so it can be used in both directions
- Height rod can be taken apart for compact storage
- Overall dimensions max. W×D×H 890×330×100 mm
- Net weight approx. 0,7 kg
- KERN MSB 80, € 80,-

5 KERN MSC 100

Portable mechanical height rod

Features

- Approved as a medical device according to 93/42/EEC
- For infants up to a maximum of 100 cm
- Measuring range 6,5–100 cm (2½–39¼ inch); increments of 1 mm (1⁄16 inch)
- Readability on scale with moveable stop
- Sturdy aluminum profile
- Stop can be folded for compact storage
- Overall dimensions max. W×D×H 1040×290×55 mm
- Net weight approx. 0,8 kg
- KERN MSC 100, € 105,-

6 KERN MSF 200

Mechanical height rod

Features

- Approved as a medical device according to 93/42/EEC
- Measuring range 60–205 cm (23½–78¾ inch); increments of 1 mm (1⁄16 inch)
- 7 Readability on scale with moveable, foldable stop or mounting on KERN scales MPS-PM (page 14), MPB-P (page 16) or 8 wall-mounted
- Sturdy aluminum profile
- Overall dimensions max. W×D×H 60×330×1630 mm
- Net weight approx. 0,7 kg
- KERN MSF 200, € 95,-

STANDARD



Elegantly-shaped baby scale – with EC type approval class III and approval for professional mobile or stationary medical use in medical diagnostics





Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Popular in paediatrics: the large tray, with its optimal ergonomic curvature, is ideal for the safe weighing of babies
- **2** Special function for neonatology centres: Difference function for calculating the weight gain before and after feeding
- Robust ABS plastic. Easy and hygienic cleaning
- Secure and non-slip positioning with rubber feet
- Hold function: When patients do not lie completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the young patient first and then take the weight reading in peace
- KERN MBC: Protective working cover included with delivery
- KERN MBC-E: Identical in construction to the MBC but without mains adapter, at a particularly attractive price, especially suitable for mobile use (eg. for midwives)

Technical data

- Backlit LCD display, digit height 25 mm
- Dimensions of tray
W×D×H 600×280×55 mm
- Overall dimensions
W×D×H 600×407×120 mm
- Overall dimensions incl. optional height rod
W×D×H 890×470×175 mm
- KERN MBC: Battery operation possible, 6×1,5 V AA, not included, operating time up to 20 h
- MBC-E: Batteries included, 6×1,5 V AA, operating time up to 20 h
- Net weight approx. 5,5 kg

Accessories

- **1** Mechanical height rod to fix to the scale, measuring range 80 cm, increments of 1 mm, KERN MBC-A01, **€ 110,-**
- **3** Foot/head rest to fix to the lying surface for even more safety when weighing small children, scope of delivery: 1 head rest and 1 foot rest, KERN MBC-A05, **€ 85,-**
- **4** Carrying case for mobile use, only for models without height rod, KERN MBC-A02, **€ 50,-**
- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- Internal rechargeable battery pack, operating time up to 40 h, without backlight, charging time approx. 12 h, only in combination with KERN YKA-44, KERN MBC-A08, **€ 54,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**


STANDARD	OPTION	FACTORY
MOVE BATT MULTI 1 DAY	ACCU MULTI MBC-E	M +3 DAYS

* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range [Max] kg	Readability [d] kg	Verification value [e] kg	Price excl. of VAT ex works €	Mandatory by law Verification	
					KERN	€
KERN						
MBC 15K2DM	6 15	0,002 0,005	0,002 0,005	285,-	965-128.	32,-
MBC 15K2DEM	6 15	0,002 0,005	0,002 0,005	265,-	965-128.	32,-
MBC 20K10M	20	0,01	0,01	285,-	965-128.	32,-
MBC 20K10EM	20	0,01	0,01	265,-	965-128.	32,-

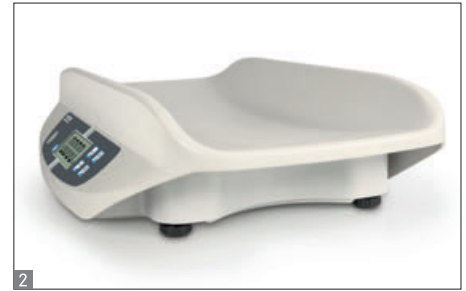Compact baby scale with special Safe Guard tray with EC type approval class III and approval for professional mobile or stationary medical use in medical diagnostics





Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Modern design
- **2** Safe-Guard tray: the special, particularly recessed tray with raised side bolsters offers optimum conditions for the safe weighing of babies
- Large, insensitive plastic tray, without corners and edges, therefore particularly hygienic and easy to clean
- Lightweight, compact design, therefore easy to transport and store
- Secure and non-slip positioning with height-adjustable rubber feet
- Easy to use
- Data hold function: for easy reading of the weight display
- Integrated AUTO-OFF function to preserve the batteries

Technical data

- LCD display, digit height 25 mm
- Dimensions of tray W×D 600×260 mm
- Overall dimensions W×D×H 600×360×141 mm
- Batteries included, 6×1,5 V AA, operating time up to 20 h
- Net weight approx. 4 kg

Accessories

- **1** Mechanical height rod to fix to the scale, measuring range 80 cm, increments of 1 mm, KERN MBA-A01, **€ 110,-**
- **1** WiFi data interface for wireless data transmission from balance to PC, tablet or smartphone, KERN YMI-A01, **€ 90,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**

STANDARD



OPTION



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range	Readability	Verification value	Price excl. of VAT ex works €	Mandatory by law Verification	
					M KERN	€
MBA 10K-3M	[Max] kg 15	[d] kg 0,005	[e] kg 0,005	425,-	965-128.	32,-



Ergonomically-optimised baby scale with large lying surface

Features

- Modern design
- Large, ergonomically optimised tray for safe weighing of neonates
- Robust plastic tray
- Easy and hygienic cleaning
- Particularly flat and lightweight construction, so it is easy to transport and store
- Secure and non-slip positioning with rubber feet
- Easy to use
- Data hold function: for easy reading of the weight display
- Special function for neonatology centres: Difference function for calculating the weight gain before and after feeding
- Internal weight memory for one weight value
- Integrated AUTO-OFF function to preserve the batteries
- Weight displayed in kg or lb
- **Note:** The models are only delivered in a set of 5 units. I.e. the price given in the table refers to a delivery of 5 items. Cannot be delivered individually.

Technical data

- LCD display, digit height 30 mm
- Dimensions of tray W×D 560×225 mm
- Overall dimensions W×D×H 650×350×55 mm
- Batteries included, 2×1,5 V AAA, operating time up to 20 h
- Net weight approx. 2 kg

Accessories

- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- **1** Carrying bag for mobile use, a carrying handle and shoulder strap ensure comfortable handling
Note: This model is only delivered in a set of 5 units. I.e. the price given refers to a delivery of 5 items. Cannot be delivered individually, KERN YTB-01S05, **€ 150,-**
- **2** Backpack for mobile use, 2 padded carrying straps and carrying handles for safe, convenient handling, additional, large zipped compartment for storing objects.
Note: This model is only delivered in a set of 5 units. I.e. the price given refers to a delivery of 5 items. Cannot be delivered individually, KERN YTB-02S05, **€ 180,-**

STANDARD



Model	Weighing range	Readability	Price set of 5 units excl. of VAT ex works €
KERN	[Max] kg	[d] kg	
MBD 20K-2S05	15	0,01	380,-

DE

KERN[®]
OPTICS

MIKROSKOPE & REFRAKTOMETER
für Labor, Industrie und die Food-Branche

PROFESSIONAL MEASURING
2021

**NEUER WEBSHOP
ONLINE!**
www.kern-sohn.com

For all details on our ranges, please see **www.kern-sohn.com**
or request your own personal catalogues right now:
Telephone +49 7433 9933-0



KERN offers you a complete, carefully-designed range of biological microscopes, stereo microscopes, metallurgical microscopes, polarisation microscopes, especially flexible, preconfigured stereomicroscope sets, individually configurable stereomicroscope systems as well as microscope cameras and camera software

Quick and easy determination of measured values, whether on site, in production or in the laboratory, ensuring safety, consistent quality and accelerating processes. This is what KERN's analog and digital refractometers stand for and at an extremely attractive price in consistently high quality.



A digital microscope with a white and black body. It features a large tablet mounted on top that displays a microscopic image of a textured surface with the word "KERN" overlaid. The microscope has two eyepieces on the left, a central objective lens, and a stage with a sample. The brand name "KERN" is visible on the side of the main body.

KERN OBE-13

KERN OBE-13

Professional personal floor scale with BMI function – with EC type approval class III and approval for professional medical use in medical diagnostics



KERN **MPE-P**
with stand

KERN **MPE-H**
with stand and
height rod



Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- **1** Practical second display on the back of the balance, which makes it easy for the patient and the doctor or nursing staff to read the weight
- **2** The MPE-HM model has an elegant height rod integrated into the column, which can be pulled out for use and has a sturdy plastic head rest and a measuring range from 88–205 cm
- The ergonomically optimised shape and robust design mean that this range is ideal for daily use in a professional environment
- Secure and non-slip positioning with height-adjustable rubber feet
- Level indicator to level the balance precisely
- Large platform with non-slip and wear-resistant surface made of plastic, to provide a secure footing for patients
- Hold function:
When patients do not stand completely still, a stable weight is calculated using an average weight and this is then “frozen”. This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- Mother and child function: the mother’s weight will be taken and then the scale will be reset to zero at the push of a button. Then the mother and baby can be weighed together. The weight of the child will now be shown on the display. This is particularly good for infants, as this way they can be cradled, safe and sound, in the arms of their mother while being weighed

- BMI function to determine underweight/normal weight/surplus weight
- High resolution readability: readability [d] for 5 sec. increased by one decimal place at the touch of a key
- **3** Easy to transport thanks to the two integrated transport wheels
- Battery- or mains-powered, rechargeable battery operation optional
- KERN MPE: Protective working cover included with delivery
- KERN MPE-E: identically with MPE though without mains adapters, only battery operation. With a particularly attractive price

Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions of weighing plate W×D×H 365×360×80 mm
- Overall dimensions
MPE-HM: W×D×H 365×570×2134 mm (with extracted height rod),
MPE-PM: W×D×H 365×570×1030 mm
- KERN MPE: Battery operation possible, 6×1,5 V AA, not included, operating time up to 20 h
- KERN MPE-E: Batteries included, 6×1,5 V AA, operating time up to 20 h
- KERN MPE: Mains adapter external, standard
- Data interface RS-232 standard
- Net weight approx. 12 kg

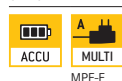
Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MPE-A01S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Internal rechargeable battery pack, operating time up to 40 h, without backlight, charging time approx. 12 h, KERN MBC-A08, **€ 54,-**
- Mains adapter external, 100–240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- Mains adapter external, 100–240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**
- Interface cable RS-232 to connect an external device, KERN CFS-A01, **€ 46,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet

STANDARD



OPTION




FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range [Max] kg	Readability [d] kg	Verification value [e] kg	Price excl. of VAT ex works €	Mandatory by law Verification	
					KERN	€
KERN						
MPE 200K-1PEM	250	0,1	0,1	400,-	965-129	53,-
MPE 250K100PM	250	0,1	0,1	430,-	965-129	53,-
MPE 200K-1HEM	250	0,1	0,1	430,-	965-129	53,-
MPE 250K100HM	250	0,1	0,1	460,-	965-129	53,-

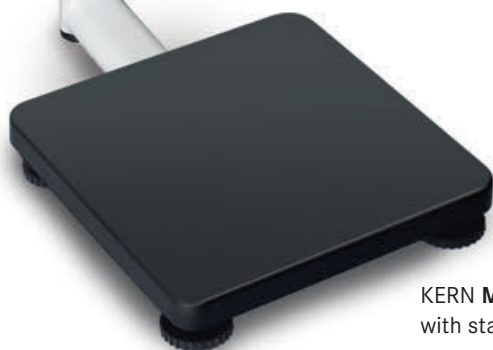
Professional personal floor scale with BMI function – with EC type approval class III and approval for professional medical use in medical diagnostics



KERN **MPS-PM**
with stand



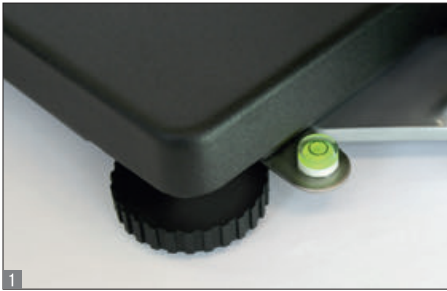KERN **MPS-M**
with free-standing
display device



Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Thanks to its compact dimensions and the good price/performance ratio, this range is ideally suited for use in treatment rooms in doctor's surgeries and hospitals
- Ergonomic display device with large numerical keypad and high-contrast LCD display for easy entry and reading of patient data, such as, for example, the height for calculating the BMI
- Sturdy steel weighing plate with a non-slip and wear-resistant surface
- Easy and hygienic cleaning
- **1** Secure and non-slip positioning with height-adjustable rubber feet
- **1** Level indicator to level the balance precisely
- Hold function:
When patients do not stand completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace

- BMI function to determine underweight/normal weight/surplus weight
- **2** Wall mount for display device, standard
- Battery- or mains-powered, rechargeable battery operation optional
- Protective working cover included with delivery

Technical data

- LCD display, digit height 25 mm
- Dimensions of weighing plate W×D×H 275×295×58 mm
- Dimensions of display device W×D×H 210×110×50 mm
- Only MPS-PM:
Overall dimensions (incl. stand) W×D×H 275×460×1015 mm
- Cable length of display device (only for models without a stand) approx. 120 cm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 50 h
- Mains adapter external, standard
- Data interface RS-232 standard
- Net weight MPS-M approx. 4.8 kg, MPS-PM approx. 6 kg

Accessories

- **3** Mechanical Height Rod, measuring range: 60–205 cm, can be screwed onto the column from the front or the back (MPS-PM) or wall-mounted, approved as a medical device according to 93/42/EEC, KERN MSF 200, **€ 95,-**
- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A02S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Internal rechargeable battery pack, operating time up to 35 h, charging time approx. 24 h, KERN FOB-A08, **€ 55,-**
- Interface cable RS-232 to connect an external device, KERN MPS-A08, **€ 46,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet

STANDARD



OPTION



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range [Max] kg	Readability [d] kg	Verification value [e] kg	Price excl. of VAT ex works €	Mandatory by law Verification	
					KERN	€
KERN						
MPS 200K100M	200	0,1	0,1	370,-	965-129	53,-
MPS 200K100PM	200	0,1	0,1	475,-	965-129	53,-



KERN MPB-P
with stand



KERN MPB
with free-standing
display device



Features

- Sturdy steel weighing plate with a non-slip and wear-resistant surface
- Easy and hygienic cleaning
- Secure and non-slip positioning with rubber feet
- Hold function:
When patients do not stand completely still, a stable weight is calculated using an average weight and this is then “frozen”. This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- Simple and convenient 4-key operation
- BMI function to determine underweight/normal weight/surplus weight
- Wall mount for the display device, only for models without a stand
- Battery- or mains-powered
- Protective working cover included with delivery

Technical data

- LCD display, digit height 25 mm
- Dimensions of weighing plate
W×D×H 315×300×60 mm
- Dimensions of display device
W×D×H 210×110×50 mm
- Only MPB-P: Overall dimensions
(incl. stand) W×D×H 315×440×1010 mm
- Cable length of display device (only for models without a stand) approx. 120 cm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 170 h
- Mains adapter external, standard
- Net weight MPB approx. 4.6 kg,
MPB-P approx. 6.4 kg

Accessories

- Mechanical Height Rod, measuring range: 60–205 cm, can be screwed onto the column from the front or the back (MPS-P) or wall-mounted, approved as a medical device according to 93/42/EEC, KERN MSF 200, **€ 95,-**
- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A02S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**

STANDARD



Model	Weighing range [Max] kg	Readability [d] kg	Price excl. of VAT ex works €
KERN			
MPB 300K100	300	0,1	215,-
MPB 300K100P	300	0,1	285,-

Professional personal floor scale – with EC type approval class III and approval for professional medical use in medical diagnostics

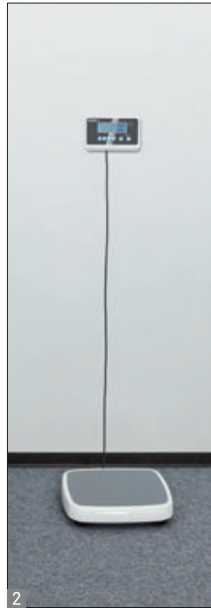
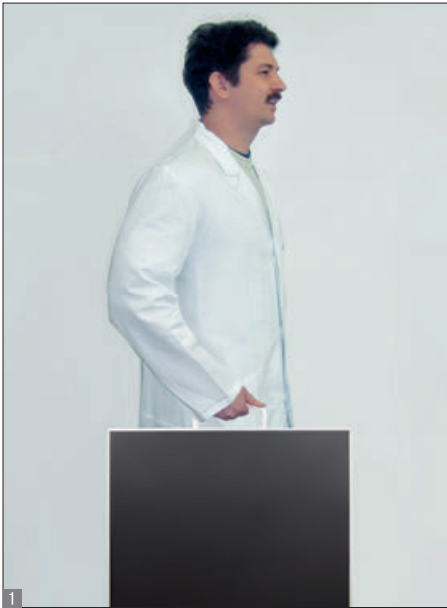


KERN MPC 250K100M

Special obesity model
MPC 300K-1M



Special obesity model
KERN MPC 300K-1LM



Features

- Version with [Max] = 300 kg and larger weighing plate! Specially suitable for weighing overweight patients who are in the obese category
- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- The bestseller for doctor's surgeries, chemists and hospitals
- Robust design for daily use in a professional environment
- Easy and hygienic cleaning
- Secure and non-slip positioning with height-adjustable rubber feet
- Level indicator to level the balance precisely
- **1** KERN MPC 300K-1LM: Sturdy handle for easy transportation of the scale
- **2** Flexible positioning of display device, e.g. free-standing or mounted to the wall. Wall mount for display device standard
- **3 + 4** Large platform with non-slip and wear-resistant surface made of plastic, to provide a secure footing for patients
- High resolution readability: readability [d] for 5 sec. increased by one decimal place at the touch of a key
- BMI function to determine underweight/normal weight/surplus weight
- Hold function:
When patients do not stand completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- Mother and child function: the mother's weight will be taken and then the scale will be reset to zero at the push of a button. Then the mother and baby can be weighed together. The weight of the child will now be shown on the display. This is particularly good for infants, as this way they can be cradled, safe and sound, in the arms of their mother while being weighed
- **5** GLP/ISO report, printout of weight value, BMI, etc. with date and time. FLP/GMP-compliant print form in DE, GB, FR, only in combination with KERN MPC-A01, CFS-A01 and printer YKB-01N
- Battery- or mains-powered, rechargeable battery operation optional
- Protective working cover included with delivery

Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions of display device
W×D×H 200×130×60 mm
- Cable length of display device approx. 150 cm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 20 h
- Mains adapter external, standard

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Internal rechargeable battery pack, operating time up to 40 h, without backlight, charging time approx. 12 h, KERN MBC-A08, **€ 54,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**
- RS-232 data interface, must be ordered at purchase, KERN MPC-A01, **€ 90,-**
- Interface cable RS-232 to connect an external device, KERN CFS-A01, **€ 46,-**

STANDARD



OPTION




FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range [Max] kg	Readability [d] kg	Verification value [e] kg	Weighing plate mm	Net weight approx. kg	Price excl. of VAT ex works €	Mandatory by law Verification	
							KERN	€
KERN MPC 250K100M	250	0,1	0,1	365×370×80 3	9	310,-	965-129	53,-
KERN MPC 300K-1M	300	0,1	0,1	365×370×80 3	9	375,-	965-129	53,-
KERN MPC 300K-1LM	300	0,1	0,1	400×500×120 4	10	525,-	965-129	53,-

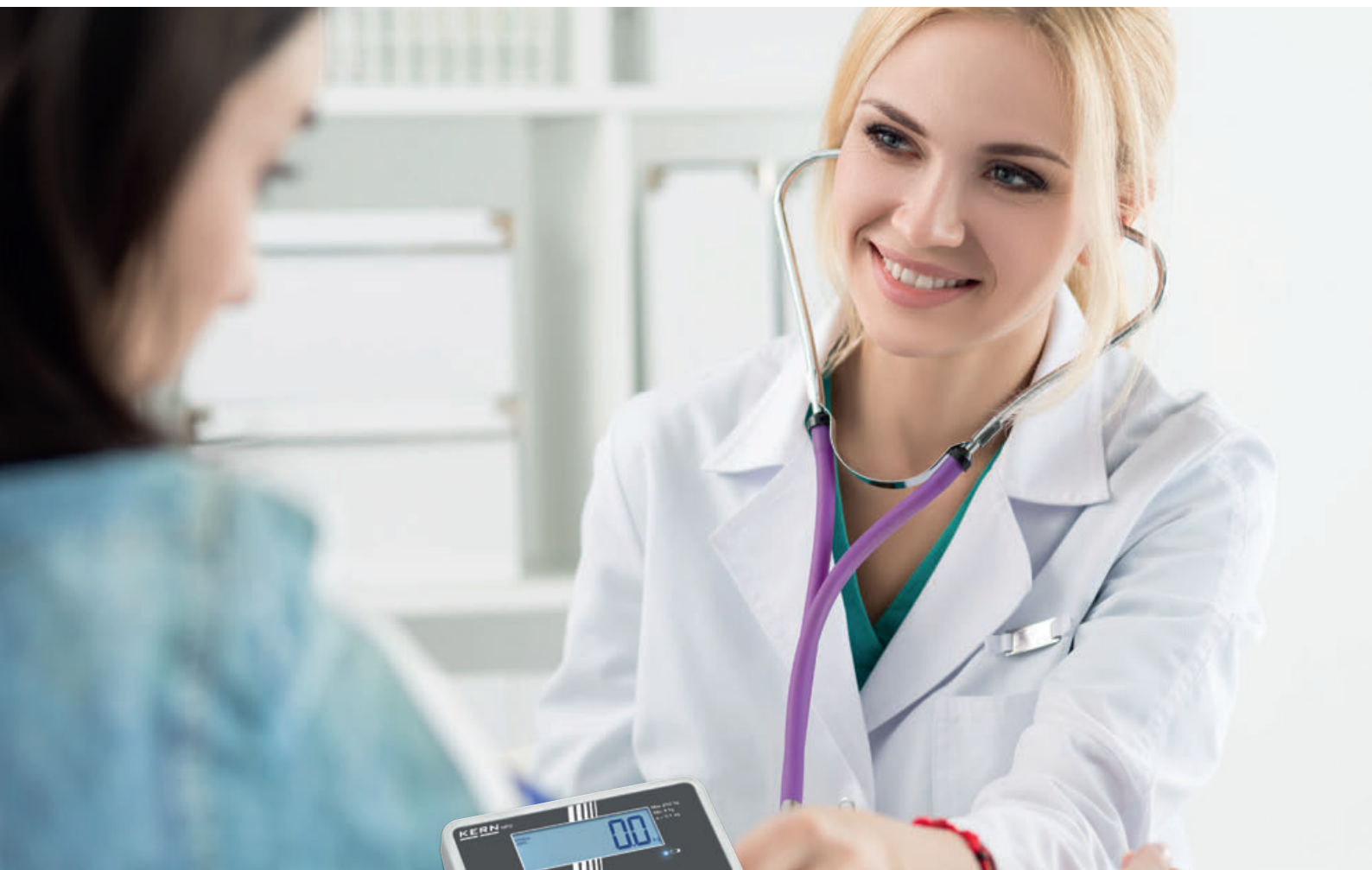Professional step-on personal floor scales approved as a medical device and with EC type approval class III



KERN MPD 250K100M



Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- **1** The bestseller for doctor's surgeries, chemists and hospitals: easy to use with no waste of time – simply tap on it, stand on it, read the weight
- Robust design for daily use in a professional environment
- Ergonomically optimised design
- Easy and hygienic cleaning
- Secure and non-slip positioning with height-adjustable rubber feet
- Level indicator to level the balance precisely
- **2** Large platform with non-slip and wear-resistant surface made of plastic, to provide a secure footing for patients
- Doctor scale MPD 200K-1EM: Identically to MPD 250K100M though without backlight and mains adapter, battery operation only, due to the reduced resolution faster display speed. Especially attractive in price!
- KERN MPD: Battery- or mains-powered, rechargeable battery operation optional
- Protective working cover included with delivery

Technical data

- Dimensions of weighing plate W×D×H 365×360×80 mm
- Dimensions of display device W×D×H 200×130×60 mm
- Net weight approx. 10 kg

MPD-E

- Large LCD display, without backlight, digit height 25 mm
- Battery operation, 6×1,5 V AA, operating time up to 20 h

MPD

- Large backlit LCD display, digit height 25 mm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 20 h
- Mains adapter external, standard

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- **3** Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- KERN MPD: Rechargeable battery pack internal, operating time up to 40 h, without backlight, charging time approx. 12 h, KERN MBC-A08, **€ 54,-**
- KERN MPD: Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- KERN MPD: Mains adapter external, 100-240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**


STANDARD	OPTION	FACTORY
BATT MULTI 1 DAY MPD	ACCU MPD	M +3 DAYS *

* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range [Max] kg	Readability [d] kg	Verification value [e] kg	Price excl. of VAT ex works €	Mandatory by law Verification	
					KERN	€
MPD 250K100M	250	0,1	0,1	290,-	965-129	53,-
MPD 200K-1EM	250	0,2	0,2	275,-	965-129	53,-



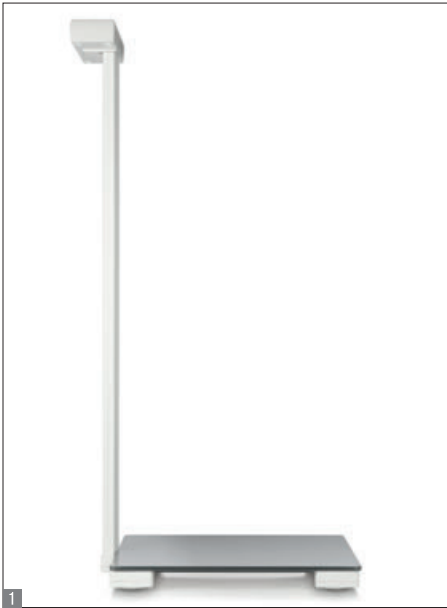KERN MPK 200K-1P



Features

- This personal scale is the ideal measuring instrument for elderly care, paramedics, fitness clubs or in the private sector and can weigh overweight persons up to 250 kg
- The column, which is firmly fixed to the scale, makes it easy to operate the scale as well as to read the weighing results even under poor lighting conditions
- Sturdy glass weighing plate for secure, comfortable standing position
- Easy and hygienic cleaning
- **1** Particularly flat design facilitates access
- Easy to operate:
simply tap – stand – read the weight
- Integrated AUTO-OFF function to preserve the batteries
- KERN MPL: BMI function to determine underweight/normal weight/surplus weight
- Weight displayed in kg or lb
- **2** Large high-contrast display that is easy to read
- KERN MPL: Adjusting function and tare function to zero the balance
- **3** Delivered in appealing packaging
- Adjusting program CAL for quick setting of the balance accuracy using an external test weight

Technical data

- LCD display, digit height 28 mm
- Dimensions of weighing plate
W×D×H 350×350×45 mm
- Overall dimensions (incl. stand)
W×D×H 350×420×955 mm
- Batteries included, 4×1,5 V AA, operating time up to 20 h
- Net weight approx. 3,8 kg

Accessories

- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**

STANDARD



Model	Weighing range [Max] kg	Readability [d] kg	Price excl. of VAT ex works €
KERN			
MPK 200K-1P	250	0,1	145,-
MPL 200K-1P	250	0,1	160,-



Personal scale with particularly large, flat weighing plate to make it easier to step on and give maximum safety

Features

- This personal scale is the ideal measuring instrument for elderly care, paramedics or privat sector, and can weigh overweight persons up to 200 kg
- Very large weighing plate
- Easy and hygienic cleaning
- **1** The scale can be easily transported using a handle and does not require much storage space
- **2** Particularly flat design facilitates access
- Very large platform with non-slip and wear-resistant plastic-surface
- Easy to operate:
simply tap – stand – read the weight
- Integrated AUTO-OFF function to preserve the batteries
- Weight displayed in kg or lb
- Large high-contrast display that is easy to read
- Delivered in appealing packaging
- **Note:** The models are only delivered in a set of 5 units. I.e. the price given in the table refers to a delivery of 5 items. Cannot be delivered individually.

Technical data

- LCD display, digit height 30 mm
- Overall dimensions
W×D×H 433×370×35 mm
- Batteries included, 4×1,5 V AAA, operating time up to 20 h
- Net weight approx. 3 kg

Accessories

- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**

STANDARD



Model	Weighing range	Readability	Price set of 5 units excl. of VAT ex works €
KERN	[Max] kg	[d] kg	
MPI 200K-1S05	200	0,1	360,-



The robust classic for quick determination of body weight

Features

- Timeless design
- Robust construction
- 1 Large, clear circular scale which is easy to read
- 2 Large platform with non-slip surface made of plastic
- Secure and non-slip positioning with rubber feet
- Easy to operate: stand – read the weight
- Easy adjustment of the display using the setting wheel
- Delivered in appealing packaging
- **Note:** The models are only delivered in a set of 5 units. I.e. the price given in the table refers to a delivery of 5 items. Cannot be delivered individually.

Technical data

- Digit height 8 mm
- Overall dimensions W×D×H 285×425×85 mm
- Net weight approx. 3,1 kg

Accessories

- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, € 35,-

STANDARD



Model	Weighing range	Readability	Price set of 5 units excl. of VAT ex works €
KERN	[Max] kg	[d] kg	
MGC 100K-1S05	150	1	290,-



Elegant look combined with large, flat glass weighing plate to make it easier to step on and give maximum safety

Features

- Elegant design
- Sturdy glass weighing plate
- MGD-L: particularly large weighing plate for secure, comfortable standing position
- Easy and hygienic cleaning
- Particularly flat design facilitates access
- Generous weighing range
- Secure and non-slip positioning with rubber feet
- Easy to operate: simply tap – stand – read the weight
- Integrated AUTO-OFF function to preserve the batteries
- Weight displayed in kg or lb
- **1** Large high-contrast display that is easy to read
- **2** Delivered in appealing packaging
- **Note:** The models are only delivered in a set of 5 units. I.e. the price given in the table refers to a delivery of 5 items. Cannot be delivered individually.

Technical data

MGD 100K-1

- LCD display, digit height 40 mm
- Batteries included, 2×CR2032, operating time up to 100 h

MGD 200K1L

- LCD display, digit height 50 mm
- Batteries included, 3×1,5 V AAA, operating time up to 100 h

Accessories

- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**

STANDARD



Model	Weighing range	Readability	Weighing plate	Net weight	Price set of 5 units excl. of VAT ex works €
KERN	[Max] kg	[d] kg	mm	kg	
MGD 100K-1S05	180	0,1	300×300×19	1,5	99,-
MGD 200K-1LS05	250	0,1	380×330×27	2	149,-



Design body analysis scale to determine body weight, body fat, body fluid, muscle mass, proportion of bone, BMI and metabolic rate/metabolic rate when exercising

Features

- Used to determine body fat, body water, muscle mass with an accuracy of 0.1 % as well as the body weight and bone mass in 0.1 kg increments
- Display showing the calorie requirements for maintaining the current body weight depending on the level of activity and sex
- Elegant design
- Sturdy glass weighing plate with four large contact surfaces for maximum precision
- Easy and hygienic cleaning
- Particularly flat design facilitates access
- Secure and non-slip positioning with rubber feet
- Easy to operate:
simply tap – stand – read the weight
- Memory for up to 10 people (gender, age, height, level of activity)
- Integrated AUTO-OFF function to preserve the batteries
- Weight displayed in kg or lb
- 1 Large high-contrast display that is easy to read
- 2 Delivered in appealing packaging
- **Note:** The models are only delivered in a set of 5 units. I.e. the price given in the table refers to a delivery of 5 items. Cannot be delivered individually.

Technical data

- LCD display, digit height 35 mm
- Overall dimensions
W×D×H 300×300×25 mm
- Batteries included, 2×CR2032, operating time up to 100 h
- Net weight approx. 2 kg

Accessories

- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, € 35,-

STANDARD



Model	Weighing range	Readability	Price set of 5 units excl. of VAT ex works €
KERN	[Max] kg	[d] kg	
MFB 150K100S05	182	0,1	140,-

Wheelchair platform scale KERN MWA

Wheelchair platform scale with two integrated access ramps for easy access – with EC type approval class III and approval for professional, stationary medical use in medical diagnostics





KERN MWA 300K-1PM



4 KERN MWA-A02: Handrail, without column, cannot be folded, for self-assembly on the KERN MWA 300K-1M wheelchair platform scale or MWA 300K-1PM. With the practical optional handrail set the KERN MWA 300K-1M wheelchair platform or MWA 300K-1PM can be converted quickly and easily into a multi-function scale. This allows for the patients to be weighed while standing on the scale, with support from the handrail or weighed in a wheelchair

Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Especially suitable for weighing patients in wheelchairs, because the low-profile platform can be approached from either side
- **1** The large steel weighing plate also makes it ideal for weighing obese patients
- Level indicator to level the balance precisely
- Hold function:
When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- **MWA 300K-1PM:** Wheelchair platform scale with practical handrail incl. column for easy, ergonomic mounting of the display device. The handrail offers patients that extra level of safety when driving or stepping onto the weighing platform
- **2** Both access ramps and the handrail KERN MWA-A01 can be folded up quickly and easily, e.g. for transportation
- BMI function to determine underweight/normal weight/surplus weight
- Battery- or mains-powered, rechargeable battery operation optional
- High resolution readability: readability [d] for 5 sec. increased by one decimal place at the touch of a key
- **3** GLP/ISO report, printout of weight value, BMI, etc. with date and time. FLP/GMP-compliant print form in DE, GB, FR, only in combination with interface cable RS-232 KERN CFS-A01 and printer YKB-01N
- Protective working cover included with delivery

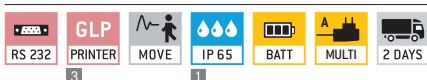
Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions of display device W×D×H 200×128×55 mm
- Overall dimensions, without handrail, W×D×H 1180×1130×90 mm
- Cable length of display device approx. 300 cm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 20 h
- Mains adapter external, standard
- Net weight approx. 76 kg

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Internal rechargeable battery pack, operating time up to 40 h, without backlight, charging time approx. 12 h, KERN MBC-A08, **€ 54,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**
- Interface cable RS-232 to connect an external device, KERN CFS-A01, **€ 46,-**
- **4** Handrail, lacquered steel, dimensions W×D×H approx. 560×40×920 mm, scope of supply: handrail, fasteners, net weight approx. 4 kg, KERN MWA-A02, **€ 85,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet

STANDARD



OPTION



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range [Max] kg	Readability [d] kg	Verification value [e] kg	Weighing surface mm	Price excl. of VAT ex works €	Mandatory by law Verification	
						KERN	€
KERN							
MWA 300K-1M	300	0,1	0,1	880×840	1140,-	965-129	53,-
MWA 300K-1PM	300	0,1	0,1	880×840	1190,-	965-129	53,-



Fig. shows
MWA 300K-1M + MWA-A04

Wheelchair scale with integrated ascending ramp and handrail set with folding patient seat – with EC type approval class III and approval for professional, stationary medical use in medical diagnostics

Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Particularly suitable for weighing patients in wheelchairs, thanks to the easy-drive-on platform with low overall height
- **1** The large steel weighing plate also makes it ideal for weighing obese patients
- Level indicator to level the balance precisely
- Hold function:
When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then “frozen”. This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- Battery- or mains-powered, rechargeable battery operation optional
- High resolution readability: readability [d] for 5 sec. increased by one decimal place at the touch of a key
- 2 handrails with folding patient seat, for self-assembly on the KERN MWA 300K-1M wheelchair platform scale, Factory Option, please order together with the balance – it can not be equipped afterwards.
- With the two practical handrails including the integrated, folding patient seat, the KERN MWA 300K-1M can easily and quickly be converted into a multi-functional scale. In this way patients can be weighed standing up, with support from the handrail or in a comfortable seated position. Thanks to an integrated damper, the seat surface can gently be lowered. When folding the patient seat up, a wheelchair can be driven onto the scale. The strut rod which is included with the delivery ensures a high level of stability and torsional rigidity. The drive-through protection for the front wheels, which is included with the delivery, gives an extra level of safety
- **2** GLP/ISO report, printout of weight value, BMI, etc. with date and time. FLP/GMP-compliant print form in DE, GB, FR, only in combination with interface cable RS-232 KERN CFS-A01 and printer YKB-01N
- Protective working cover included with delivery

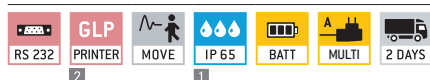
Technical data

- Dimensions W×D×H approx.
1132×880×960 mm (mounted)
- Dimensions of seating surface W×D 470×310 mm
- Height of the seat above the platform 460 mm
- Scope of supply: 2 handrails, strut rod, drive-through protection, folding seat, fasteners, Factory option, the assembly takes place in the factory KERN, the delivery is completely assembled on pallet
- Net weight approx. 84 kg

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- Internal rechargeable battery pack, operating time up to 40 h, without backlight, charging time approx. 12 h, KERN MBC-A08, **€ 54,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**
- Interface cable RS-232 to connect an external device, KERN CFS-A01, **€ 46,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet

STANDARD



OPTION



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range	Readability	Verification value	Weighing surface	Price excl. of VAT ex works €	Mandatory by law Verification	
KERN	[Max] kg	[d] kg	[e] kg	mm		M III KERN	€
MWA 300K-1M	300	0,1	0,1	880×840	1140,-	965-129	53,-
MWA-A04	Handrail set with folding patient seat, please order at the same time				770,-	-	-



Compact wheelchair platform scale with two integrated access ramps for easy access

Features

- Especially suitable for weighing patients in wheelchairs because of the low-profile platform which can be approached from either side
- **1** The large steel weighing plate also makes it ideal for weighing obese patients
- Ideal also for rapid weighing of, e.g. laundry carts, container trolleys, roller containers, etc.
- Secure and non-slip positioning with height-adjustable rubber feet
- **2** Level indicator to level the balance precisely
- Hold function:
When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then “frozen”. This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- **3** The scale can be easily transported using the two handles and two rollers and does not require much storage space
- Wall mount for the display device
- Protective working cover included with delivery

Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions of display device W×D×H 200×130×60 mm
- Overall dimensions, W×D×H 1050×1050×100 mm
- Batteries included, 6×1,5 V AA, operating time up to 40 h
- Cable length of display device approx. 300 cm
- Net weight approx. 55 kg

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- For more details see the internet

STANDARD



Model	Weighing range	Readability	Weighing surface	Price excl. of VAT ex works €
KERN	[Max] kg	[d] kg	mm	
MWB 300K-1	300	0,1	800×800	690,-

Wheelchair platform scale and stretcher scale with low overall height for easy access – with EC type approval class III and approval for professional, mobile medical use in medical diagnostics



1 KERN MWS



2 KERN MWS-L



Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Especially suitable for weighing patients in wheelchairs because of the low-profile platform which can be approached from either side and the memory function for tare weights for up to five (power) wheelchairs
- **1** KERN MWS: Three balances in one:
With the practical KERN MWS-A02 handrail set, the wheelchair platform scale MWS 300K100M can be converted quickly and easily into a multi-function scale. In this way patients can be weighed standing up on their own, using the support of the handrail or in a wheelchair. For a more ergonomic system the display device can be fitted directly onto the handrail
- **2** KERN MWS-L: Due to its large dimensions and large weighing range, it is especially suitable for weighing immobile patients on stretchers or wheelchairs
- The large steel weighing plate also makes it ideal for weighing obese patients
- Secure and non-slip positioning with height-adjustable rubber feet
- Level indicator to level the balance precisely

- Hold function: When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- Can easily be transported using 2 rollers and a handle and also does not require much storage space.
- Display device with two magnets on the back, ideal to fix it to metallic surfaces, e.g. to the weighing plate during transportation
- Battery- or mains-powered, rechargeable battery operation optional
- Protective working cover included with delivery

Technical data

- LCD display, digit height 25 mm
- Dimensions of display device W×D×H 210×110×50 mm
- Overall dimensions
MWS 300K100M B×T×H 1155×800×66 mm
MWS 400K100DM B×T×H 1255×1060×69 mm
MWS 300K-1LM B×T×H 1500×860×68 mm
- Overall dimensions incl. MWS-A02 W×D×H 1155×820×1070 mm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 50 h

- Mains adapter external, standard
- Cable length of display device approx. 180 cm

Accessories

- **1** Handrail set to fit to the wheelchair platform scale KERN MWS 300K100M, W×D×H approx. 1150×60×900 mm. Scope of delivery: two handrails, a stability bar, fasteners. KERN MWS-A02, **€ 550,-**
- Stand to elevate display device, height of stand approx. 950 mm, KERN MWS-A01, **€ 180,-**
- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A02S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Internal rechargeable battery pack, operating time up to 50 h, charging time approx. 24 h, KERN FOB-A08, **€ 55,-**
- Interface cable RS-232 to connect an external device, KERN MPS-A08, **€ 46,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet


STANDARD	OPTION	FACTORY
RS 232 MOVE BATT MULTI 2 DAYS	ACCU	M +3 DAYS

* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range	Readability	Verification value	Weighing surface	Net weight	Price excl. of VAT ex works €	Mandatory by law Verification	
KERN	[Max] kg	[d] kg	[e] kg	mm	approx. kg		KERN	€
MWS 300K100M	300	0,1	0,1	900×740	30	1195,-	965-129	53,-
MWS 400K100DM	300 400	0,1 0,2	0,1 0,2	1000×1000	44	1390,-	965-130	83,-
MWS 300K-1LM	300	0,1	0,1	1200×800	42	1550,-	965-129	53,-

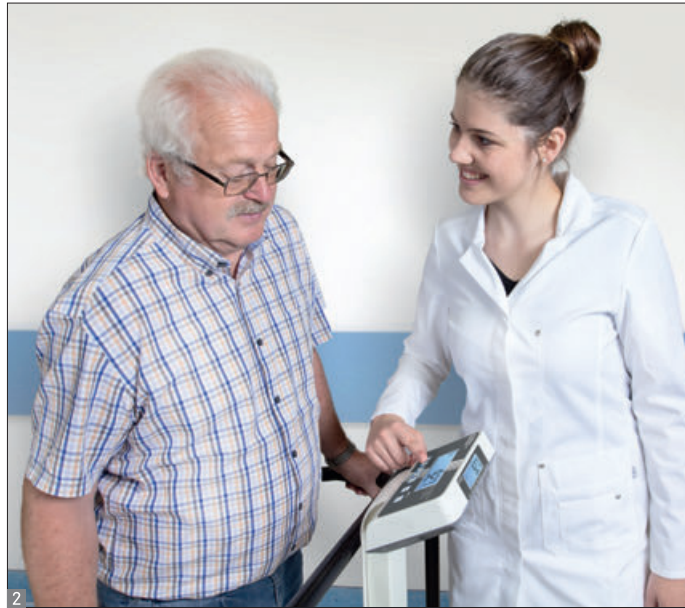Sturdy handrail scale for a secure feeling during weighing – with EC type approval class III and approval for professional medical use in medical diagnostics





Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- Ideal for care homes, centres for geriatric medicine, internal medicine and many more
- **1** The sturdy support structure helps fragile and weak patients to stand safely and, at the same time, allows a precise determination of the patients' weight
- Large surface, for the patient to have the most secure footing, hygienic and easy to clean
- Particularly easy access to the platform through its low height of just 30 mm
- Level indicator to level the balance precisely
- Hold function:
When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- **2** Practical second display on the back of the balance, which makes it easy for both the patient and the doctor or nursing staff to read the weight

- BMI function to determine underweight/normal weight/surplus weight
- Ergonomic display device with large keys and a high-contrast LCD display for easy entry and reading of patient data or the height for calculating the BMI
- Battery- or mains-powered, rechargeable battery operation optional
- Protective working cover included with delivery

Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions of weighing plate
W×D×H 780×680×68 mm
- Dimensions of weighing surface
W×D 600×600 mm
- Dimensions of display device
W×D×H 200×128×55 mm
- Overall dimensions
W×D×H 780×680×1100 mm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 50 h
- Mains adapter external, standard
- Net weight approx. 40 kg

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MPE-A01S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Rechargeable battery pack internal, operating time up to 40 h without backlight, charging time approx. 12 h, KERN MBC-A08, **€ 54,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, CH, UK, KERN YKA-43, **€ 75,-**
- Mains adapter external, 100-240 V, 50/60 Hz, standard EU, KERN YKA-44, **€ 45,-**
- Interface cable RS-232 to connect an external device, KERN CFS-A01, **€ 46,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- Further details, see Internet

STANDARD



OPTION



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range	Readability	Verification value	Price excl. of VAT ex works €	Mandatory by law Verification	
					KERN	€
KERN MTA 400K-1M	[Max] kg 300 400	[d] kg 0,1 0,2	[e] kg 0,1 0,2	1250,-	965-130	83,-

Mobile chair scale with four wheels – with EC type approval class III and approval for professional medical use in medical diagnostics





Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- **1** Mobile version with locking brakes on all wheels
- With its four wheels, this chair scale provides maximum mobility in bringing the scale to the patient. This ensures more efficient use of time for the clinic staff and greater safety for the patients who can be weighed in their familiar environment
- The four steerable rollers, with their large diameter, make it easier to get over door thresholds, edges and across the gaps into elevators
- This chair scale is the ideal measuring instrument for retirement homes, rehabilitation centres and clinics with lots of patients, and can weigh overweight persons up to 300 kg
- For fragile patients, the comfortable, ergonomically optimised seat offers secure support during weighing
- **2** Two foldable armrests and footrests make transfers into the chair easier. Ideal for overweight patients or barrier-free use, e.g. for transfers from a bed to the scale

- Hold function:
When patients do not sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- High resolution readability: readability [d] for 5 sec. increased by one decimal place at the touch of a key
- Protective working cover included with delivery

Technical data

- **3** LCD display, digit height 25 mm
- Dimensions of display device W×D×H 210×110×50 mm
- Dimensions of seat surface W×D 380×360 mm
- Overall dimensions W×D×H 647×860×910 mm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 70 h
- Mains adapter external, standard
- Net weight approx. 22 kg

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A02S05, **€ 40,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**
- Internal rechargeable battery pack, operating time up to 50 h, charging time approx. 24 h, KERN FOB-A08, **€ 55,-**
- Interface cable RS-232 to connect an external device, KERN MPS-A08, **€ 46,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet

STANDARD



OPTION



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range	Readability	Verification value	Price excl. of VAT ex works €	Mandatory by law Verification	
					KERN	€
KERN	[Max] kg	[d] kg	[e] kg			
MCB 300K100M	300	0,1	0,1	995,-	965-129	53,-

Ergonomically optimised chair scale – with EC type approval class III and approval for professional medical use in medical diagnostics





Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- This chair scale is the ideal measuring instrument for retirement homes, rehabilitation centres and clinics with lots of patients, and can weigh overweight persons up to 250 kg
- **1** Mobile version with 2 steerable rollers and particularly convenient locking brakes at the rear
- With its four wheels, this chair scale provides maximum mobility in bringing the scale to the patient. This ensures more efficient use of time for the clinic staff and greater safety for the patients, who can be weighed in their familiar environment
- The rollers, with their large diameter, make it easier to get over door thresholds, edges and across the gaps into elevators
- For fragile patients, the comfortable, ergonomically optimised seat offers secure support during weighing
- Because of the clearly contrasting black colour of the seating surface, footrests and armrests ideally suitable for dementia patients

- **2** Two foldable armrests and footrests make transfers into the chair easier. Ideal for overweight patients or barrier-free use, e.g. for transfers from a bed to the scale
- Ergonomically positioned carrying handles
- Hold function:
When patients do not sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- High resolution readability: readability [d] for 5 sec. increased by one decimal place at the touch of a key
- Protective working cover included with delivery

- Rechargeable battery pack internal, standard, operating time up to 40 h without backlight, charging time approx. 12 h
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 20 h
- **4** Mains adapter, external, included, with strain relief and disconnector plug to protect the power supply components
- Net weight approx. 24 kg

Accessories

- **5** Practical mains adapter pouch to store the mains adapter which is supplied as standard. Fitting option using two Velcro fasteners, KERN MCC-A01, € 30,-
- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, € 40,-
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, € 35,-

Technical data

- **3** Large backlit LCD display with digit height 25 mm
- Dimensions of display device
W×D×H 200×130×60 mm
- Dimensions of seat surface
W×D 360×390 mm
- Overall dimensions
W×D×H 625×990×985 mm


STANDARD



FACTORY



* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.

Model	Weighing range	Readability	Verification value	Price excl. of VAT ex works €	Mandatory by law Verification	
					 KERN	€
MCC 250K100M	[Max] kg 250	[d] kg 0,1	[e] kg 0,1	835,-	965-129	53,-

Mobile chair scale in ergonomically-optimised design for safe, convenient weighing





Features

- This chair scale is the ideal measuring instrument for elderly care and paramedics and is also suitable for overweight persons with a weight up to 300 kg
- **1** Mobile version with 4 rollers and particularly convenient locking brakes on both rear wheels
- With its four wheels, this chair scale provides maximum mobility by bringing the scale to the patient. This ensures more efficient use of time for the staff and greater safety for the patients who can be weighed in their familiar environment
- The rollers, with their large diameter, make it easier to get over door thresholds, edges and across the gaps into elevators
- For fragile patients, the comfortable, ergonomically optimised seat offers secure support during weighing
- Because of the clearly contrasting black colour of the seating surface and footrests ideally suitable for dementia patients
- **2** Two foldable armrests and footrests make transfers into the chair easier. Ideal for overweight patients or barrier-free use, e.g. for transfers from a bed to the scale
- Ergonomically positioned carrying handles
- Hold function:
When patients do not sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- Practical battery operation using standard batteries ensures a high level of flexibility and freedom from mains adapters, sockets, chargers, etc.
- High resolution readability: readability [d] for 5 sec increased by one decimal place at the touch of a key
- Protective working cover included with delivery

Technical data

- **3** LCD display, digit height 25 mm
- Dimensions of display device W×D×H 200×130×60 mm
- Dimensions of seat surface W×D 360×390 mm
- Overall dimensions W×D×H 625×990×985 mm
- Batteries included, 6×1,5 V AA, operating time up to 20 h
- Net weight approx. 20 kg

Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05, **€ 40,-**
- Practical pocket to fix on the back of the chair, e.g. for the user instructions, cleaning cloths, etc. KERN MCC-A01, **€ 30,-**
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01, **€ 35,-**

STANDARD



Model	Weighing range	Readability	Price excl. of VAT ex works €
KERN	[Max] kg	[d] kg	
MCD 300K-1	300	0,1	750,-



Hospital need

KERN scales and measuring instruments for clinics, medical practices and rehabilitation facilities

Reliable scales are not only required in direct contact with patients, but also in adjacent areas such as laboratories, research, pathology, hospital laundries, etc. a wide variety of scales are used. These scales also have to meet high requirements, such as robust, dirt-repellent stainless steel surfaces that are hygienic and easy to clean. Simple operation and clearly legible displays speed up daily work, save costs and help to avoid mistakes.

The practical dynamometers, which we offer in various designs and measuring ranges, are particularly suitable for measuring muscle forces in the context of rehabilitation measures.



Page 43



Stainless steel drive-through scale with two integrated access ramps and EC type approval [M]

Features

- Robust stainless steel drive-through scale for rapid weighing of e.g. laundry carts, container trolleys, roller containers, etc. Ideal for hospital laundry services, goods inwards, hospital kitchens, etc.
- Low platform height and integrated access ramps on both sides facilitate access. No need for pit frame installation – which saves money
- Weighing bridge stainless steel, extremely resistant to bending
- **1** 4 load cells steel, encapsulated. Protection against dust and water splashes IP68, suitable for continuous use in wet areas
- **2** Stainless steel display device with protection against dust and water splashes IP65
- Superior display size: digit height 52 mm. Bright backlight for easy reading of weight, even in poor lighting conditions
- Easy and hygienic cleaning
- Suitable for the ever-increasing hygienic requirements in the medical environment
- Totalising of weights and piece counts

Technical data

- Large backlit LCD display, digit height 52 mm
- Overall dimensions
W×D×H 1600×1200×80 mm
- Weighing surface, without access ramps
W×D 1000×1000 mm
- Dimensions of display device
W×D×H 266×165×96 mm
- Cable length of display device approx. 5 m
- Permissible ambient temperature
-10 °C/40 °C

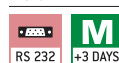
Accessories

- Rechargeable battery pack internal, operating time up to 35 h without backlight, charging time approx. 10 h, must be ordered at purchase, KERN GAB-A04, **€ 40,-**
- **3** Stand to elevate display device, stainless steel, height-adjustable, height of stand 750–1000 mm, KERN BFN-A04, **€ 285,-**
- Pair of base plates to fix the weighing bridge to the floor, KERN BFN-A03, **€ 60,-**
- Large display with superior display size, only in combination with KFN-A01, please refer to KERN's website for more details. KERN YKD-A02, **€ 220,-**
- Cable with special length 15 m, between display device and platform, must be ordered at purchase, KERN BFB-A03, **€ 160,-**
- Data interface RS-232, interface cable included, approx. 1.5 m, must be ordered at purchase, KERN KFN-A01, **€ 95,-**
- Matrix needle printer, KERN YKN-01, **€ 230,-**
- Label printer, KERN YKE-01, **€ 590,-**
- Thermal printer, KERN YKB-01N, **€ 310,-**
- For further details, plenty of further accessories and suitable printers see Internet

STANDARD



FACTORY



Model	Weighing range	Readability	Verification value	Minimum load	Net weight	Price excl. of VAT ex works €	Mandatory by law Verification	
KERN	[Max] kg	[d] g	[e] g	[Min] g	kg		KERN	€
NFN 600K-1M	600	200	200	4000	125	2650,- ↓	965-230	132,-
NFN 1.5T-4M	1500	500	500	10000	135	2850,- ↓	965-230	132,-

↓ Price reduction



Organ scale stainless steel and IP68 protected – now with improved accuracy

Features

- Particularly suitable to determine the specific weight of organs, tissues, etc., in pathology or in the laboratory
- Tough industry standard suitable for use in harsh industrial applications
- IP68 protection: very high type of protection against dust and water. Making it very robust and durable. Ideal for harsh ambient conditions. Submersion in water even possible
- Stainless steel housing, rust-proof and easy to clean
- Very fast display: steady weight values within 2 sec
- High mobility: Battery operation, compact construction and low weight make this balance ideal for use in several locations
- Carrying grips on the underside of the scale for easy transportation

Technical data

- Large LCD display, digit height 25 mm
- Dimensions of weighing plate, stainless steel W×D 236×195 mm
- Overall dimensions W×D×H 240×280×120 mm
- Ready for use: Batteries included, 4×1,5 V size D, operating time without backlight FXN-N up to 250 h, FXN-M up to 625 h, AUTO-OFF function to preserve the batteries, can be switched off
- Net weight approx. 3,2 kg
- Permissible ambient temperature -10 °C/40 °C

Accessories

- Tare pan made of stainless steel, ideal for weighing organs, tissues, etc. W×D×H 370×240×20 mm, KERN RFS-A02, **€ 55,-**

STANDARD



FACTORY



Model	Weighing range [Max] kg	Readability [d] g	Verification value [e] g	Minimum load [Min] g	Price excl. of VAT ex works €	Option Verification	
						KERN	€
KERN							
FXN 3K-4N	3	0,5	–	–	305,-	–	–
FXN 6K-3N	6	1	–	–	305,-	–	–
FXN 10K-3N	15	2	–	–	305,-	–	–
FXN 30K-3N	30	5	–	–	305,-	–	–
Models with EC type approval [M]							
FXN 3K-3M	3	1	1	20	370,-	965-227	55,-
FXN 6K-3M	6	2	2	40	370,-	965-228	70,-
FXN 10K-3M	15	5	5	100	370,-	965-228	70,-
FXN 30K-2M	30	10	10	200	370,-	965-228	70,-



Hand grip dynamometer, e.g. for rehabilitation treatment after accidents

Features

- Especially suitable for use in rehabilitation clinics for determining manual clamping force
- There are four measuring methods, for example, as part of a rehabilitation program, help the medical staff to monitor the strength of the patient's hands and carry out controlled training:
 - Real time mode: immediately shows the client's current strength
 - Peak/Max mode: shows the maximum strength of a client's grip
 - Average mode: Calculates the average strength from two grips
 - Counting mode: Counts the number of presses which exceed a previously defined strength limit
- An ideal device to determine reduced handstrength and a possible mortality risk of elderly persons as well as a malnutrition in case of chemotherapy or similar treatments.
- Safe, comfortable use thanks to non-slip rubber grips
- Integrated AUTO-OFF function after 1 minute to preserve the batteries
- Result displayed in kg or lb
- MAP 80K1S: Special version for children: The small handle depth allows children in particular to easily and ergonomically grip the handles
- MAP 130K1: Special version for body builders: Its design and extended measuring range mean that it offers additional capacity, which can accommodate the higher fundamental force exerted by body builders
- 1 Exchangeable springs facilitate fast switching of the capacity (additional spring sets are included with delivery). The varying rigidity of the individual springs makes the hand grip dynamometer ideal for a wide variety of patient groups, e.g. children or senior citizens or in sports medicine
- 2 Stable case for safe, easy transport and for storage of the additional spring sets, standard, W×D×H 350×265×85 mm

Technical data

- LCD display, digit height 12 mm
- Batteries included, 1×CR2450, operating time up to 53 h
- Net weight approx. 0.3 kg

STANDARD

UNIT


BATT

1 DAY

Model	Measuring range	Readability	Spring sets	Overall dimensions W×D×H	Price excl. of VAT ex works €	Option	
						ISO Calibr. Certificate	
KERN	[Max] kg	[d] kg	kg	mm		ISO KERN	€
MAP 80K1S	80	0,1	10, 20, 40, 80	55×88×212	235,-	961-167	120,-
MAP 80K1	80	0,1	20, 40, 80	55×102×212	235,-	961-167	120,-
MAP 130K1	130	0,1	40, 80, 130	55×102×212	235,-	961-167	120,-



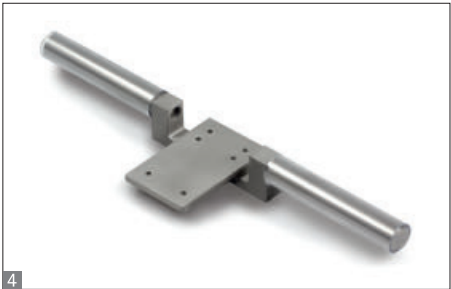See images for dynamometers with optional accessories AC 45-47



Force measurement for medical applications

Features

- Determining muscle function and force is used within many areas of medical diagnostics. Amongst others, in
 - orthopaedics for determining the function of the musculoskeletal system
 - physiotherapy for restoring motivity
 - occupational therapy for treating physical damage
- Within sport, measuring muscular force development can be used alongside training
- These measurements can be supported by the SAUTER FK (alternatively SAUTER FH) force gauge when connected with the three physio sensors (AC 45, AC 46, AC 47) as shown
- These sensors to measure muscular strength have been designed in an ergonomic manner. The surface is made of soft foam rubber, which sits comfortably on the skin
- Connecting the precision measuring device with the optional stainless steel handles SAUTER AFK 02 | AFH 04 means that you can obtain reliable weighing results
- Our recommendation: Combine the force measuring device (Sauter FK or FH) with the stainless steel handles (Sauter AFK 02 or AFH 04, depending on the force measuring device selected) and one or more physio sensors (AC 45, AC 46, AC 47)
- All components can be ordered individually

STANDARD FK



Model	Measuring range [Max] N	Readability [d] N	Price excl. of VAT ex works €
SAUTER FK 50.	50	0,02	250,-
FK 100.	100	0,05	250,-
FK 250.	250	0,1	250,-
FK 500.	500	0,2	250,-
FK 1K.	1000	0,5	250,-
FH 50.	50	0,01	460,-
FH 100.	100	0,05	460,-
FH 200.	200	0,1	460,-
FH 500.	500	0,1	460,-

STANDARD FH




Model	Description	Price excl. of VAT ex works €
SAUTER 4 AFK 02	Stainless steel handle with plastic coating for force measuring devices SAUTER FK	95,-
4 AFH 04	Stainless steel handle with plastic coating for force measuring devices SAUTER FH	95,-


OPTION FH



Model	Description	Price excl. of VAT ex works €
SAUTER 1 AC 45	Concave force sensor with optimised radius for measurement, particularly of arms and legs, up to 1 kN	135,-
2 AC 46	Flat square-shaped sensor for lateral power sensing e.g. of back, chest or arm, up to 1 kN	90,-
3 AC 47	Round sensor to measure particular muscle groups, such as, for example, the shoulder, up to 1 kN	95,-

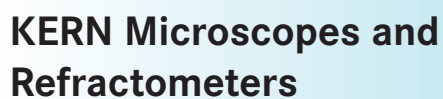

Furthermore, KERN offers a wide range of analogue and digital refractometers. Quickly and easily determine reliable measured values, no matter whether in the laboratory or on site, this ensures safety, consistent quality and accelerates the processes. This is what the analogue and digital refractometers from KERN stand for, and at an extremely attractive price in consistently high quality. The refractometers are mainly used in hospitals, doctors' surgeries, medical training facilities, nursing homes, sports medicine (doping control), veterinary surgeries, etc.

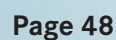



! Illustration shows optional accessories, tablet with integrated camera KERN ODC 241, further details see Internet



Elegant, dynamic and impressive – the new all-round compound microscope for schools, training and laboratories

Features

- The OBE-12/13 range stands out through its exclusive, dynamic device, which is second to none in terms of sturdy construction and ergonomics. The clever storage compartment on the back enables quick practical storage for your power cable. Thanks to the USB connection technology, it is also possible to supply power using an external powerbank
- The impressive, infinitely dimmable 3 W LED guarantees bright illumination of your sample
- A further highlight is the Butterfly tube which is integrated as standard and which enables you to achieve the ideal viewing angle. The height-adjustable and thereby focusable 1.25 Abbe condenser with aperture diaphragm is a further quality feature of the OBE range and guarantees the very best concentration of light
- Height adjustment of the fully-equipped mechanical stage is carried out using a coarse and fine focusing knob on both sides. The ergonomically designed coaxial drive enables you to work with the samples and move them rapidly
- A large selection of different eyepieces and objectives are available to you as accessories
- 1 Monocular version, 2 Trinocular version, 3 Butterfly tube
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the model outfit list, see Internet

Scope of application

- Training, haematology, sediment investigation, doctor's practise

Applications/Samples

- Translucent, thin, high-contrast, less complex samples (e.g. plant tissue, coloured cells/parasites)

Technical data

- Finite optical system
- Quadplex nosepiece
- Butterfly 30° inclined
- Diopter adjustment: One-sided (for binocular and trinocular models)
- Overall dimensions W×D×H 360×150×320 mm
- Net weight approx. 4,6 kg

STANDARD



OPTION



Model	Standard configuration					Price excl. of VAT ex works €
	Tube	Eyepiece	Objective quality	Objectives	Illumination	
KERN						
OBE 121	Monocular	HWF 10×/ø 18 mm	Achromatic	4×/10×/40×	3W LED (transmitted)	425,-
OBE 122	Binocular	HWF 10×/ø 18 mm	Achromatic		3W LED (transmitted)	520,-
OBE 124	Trinocular	HWF 10×/ø 18 mm	Achromatic		3W LED (transmitted)	630,-
OBE 131	Monocular	HWF 10×/ø 18 mm	Achromatic	4×/10×/40×/100×	3W LED (transmitted)	490,-
OBE 132	Binocular	HWF 10×/ø 18 mm	Achromatic		3W LED (transmitted)	580,-
OBE 134	Trinocular	HWF 10×/ø 18 mm	Achromatic		3W LED (transmitted)	690,-



Also available as digital, trinocular and phase contrast model. KERN OBL-1 also available with EPI fluorescence illumination unit on request



The high-performance compound microscope for every laboratory, hospital and doctor's practice with fixed, pre-centred Koehler illumination

Features

- The KERN OBF-1 and OBL-1 models are excellent, stable laboratory microscopes for all common routine applications
- Depending on the application, there is a choice of models with strong, continuously dimmable 3 W LED or 20 W halogen illumination (Philips)
- **1** Trinocular models as well as specially pre-configured phase contrast models available, please feel free to contact us
- The fixed, pre-centred and focusable 1,25 Abbe condenser with aperture diaphragm and field diaphragm gives you a simplified Koehler illumination, the centre can not be moved
- The large mechanical stage and its specimen holder holds up to two samples at the same time, which can be focused quickly and easily by a coaxial coarse and fine drive on both sides
- A large selection of eyepieces, **2** objectives and colour filters as well as a darkfield condenser, a simple **3** polarising unit, different **4** phase contrast kits through to HBO and LED fluorescence units are available to you as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the model outfit list

Please refer to the website for further details as well as information on the full range of model features

Scope of application

- Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, oncology, entomology, vets, water analysis, breweries, sewage plants

Applications/Samples

- Translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissues)

Technical data

- Finite optical system DIN (OBF)/Infinity optical system (OBL)
- Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- Diopter adjustment: One-sided
- Overall dimensions W×D×H 395×200×380 mm
- Net weight approx. 6,7 kg

STANDARD



OPTION



Model	Standard configuration					Price excl. of VAT ex works €
	Tube	Eyepiece	Objective quality	Objectives	Illumination	
KERN						
OBF 121*	Binocular	HWF 10×/ø 18 mm	Achromatic	4×/10×/ 40×/100×	6 V/20 W Halogen (transmitted)	720,- ↓
OBF 122*	Binocular	HWF 10×/ø 18 mm	Plan		6 V/20 W Halogen (transmitted)	830,- ↓
OBF 123*	Binocular	HWF 10×/ø 18 mm	Plan		3 W LED (transmitted)	810,- ↓
OBF 131*	Trinocular	HWF 10×/ø 18 mm	Achromatic		6 V/20 W Halogen (transmitted)	830,- ↓
OBF 132*	Trinocular	HWF 10×/ø 18 mm	Plan		6 V/20 W Halogen (transmitted)	940,- ↓
OBF 133*	Trinocular	HWF 10×/ø 18 mm	Plan		3 W LED (transmitted)	920,- ↓
OBL 125	Binocular	HWF 10×/ø 20 mm	Infinity E-Plan		6 V/20 W Halogen (transmitted)	910,-
OBL 127	Binocular	HWF 10×/ø 20 mm	Infinity E-Plan		3 W LED (transmitted)	910,-
OBL 135	Trinocular	HWF 10×/ø 20 mm	Infinity E-Plan		6 V/20 W Halogen (transmitted)	1020,-
OBL 137	Trinocular	HWF 10×/ø 20 mm	Infinity E-Plan		3 W LED (transmitted)	1020,-

* ONLY WHILE STOCKS LAST!

↓ Price reduction



OBN-13



OBN-15

Professionalism and versatility united in one microscope – with Koehler illumination for demanding applications



! Also available as fluorescence microscope

Features

- The OBN series stands out because of its unbeatable and consistently high quality and its ergonomic design. The range of modular components means that the OBN series can be individually customised for the professional user
- Depending on the application, there is a choice of models with strong, continuously dimmable 3 W LED or 20 W halogen transmitted illumination (Philips)
- In addition the halogen variant is available as a pre-configured phase contrast microscope, which, through the combination of a 2 professional quintuple condenser wheel, phase contrast condenser and Infinity Plan phase contrast objectives makes it a high-quality, fully-equipped microscope for all applications related to contrast procedures
- 1 This series has a professional Koehler illumination unit with an adjustable field diaphragm as well as a height-adjustable 1,25 Abbe condenser which can be centred and which has an adjustable aperture diaphragm
- The extremely large mechanical stage with ergonomic, coaxial coarse and fine focusing knob on both sides enables you to adjust and focus your sample rapidly and accurately
- A wide variety of modular systems, such as, for example, a swing-out condenser, various eyepieces, objectives, colour filters, phase contrast units, a darkfield condenser, a simple polarising unit, Butterfly Tube, through to complete fluorescence units are available to you as accessories
- This centring eyepiece for adjusting the phase contrast (OBN 158), a protective dust cover, eye cups as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera. You can select this adapter from the model outfit list

Scope of application

- Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, sewage treatment plants, oncology, entomology, vets, water analysis and breweries

Applications/Samples

- Translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissue)

Technical data

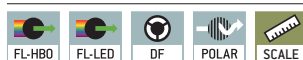
- Infinity optical system
- Quintuple nosepiece
- Siedentopf 30° inclined/360° rotatable
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H
390×200×400 mm
- Net weight approx. 9 kg

Please refer to the website for further details as well as information on the full range of model features


STANDARD



OPTION



Model	Standard configuration					Price excl. of VAT ex works €
	Tube	Eyepiece	Objective quality	Objectives	Illumination	
KERN						
OBN 132	Trinocular	HWF 10×/ø 20 mm	Infinity Plan	4×/10×/20×/ 40×/100×	6 V/20 W Halogen (transmitted)	1370,-
OBN 135	Trinocular	HWF 10×/ø 20 mm	Infinity Plan		3 W LED (transmitted)	1390,-
OBN 158	Trinocular	HWF 10×/ø 20 mm	Infinity Plan	4×/PH10×/PH20×/ PH40×/PH100×	6 V/20 W Halogen (transmitted)	2540,-
OBN 141	Trinocular	WF 10×/ø 20 mm	Infinity Plan		LED + 3 W LED Epi Fluorescence (B/G)	4470,-
OBN 147	Trinocular	WF 10×/ø 20 mm	Infinity Plan	4×/10×/20×/ 40×/100×	Halogen + 100 W Epi Fluorescence (B/G)	4390,-
OBN 148	Trinocular	WF 10×/ø 20 mm	Infinity Plan		Halogen + 100 W Epi Fluorescence (B/G/UV/V)	5630,-

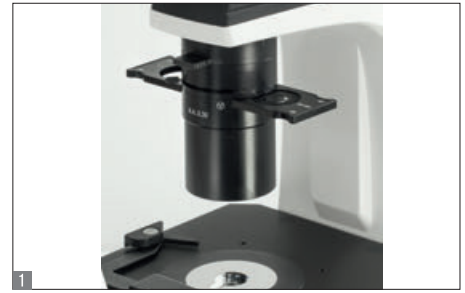

The inverted biological laboratory microscope – also with fluorescence

Features

- The OCM range stands out through its design which is ergonomic, robust and extremely stable. This design, with its large working distance, is particularly suitable for the monitoring and analysis of cell cultures, for example
- Strong and infinitely adjustable 30 W halogen illumination unit ensures the optimum illumination in the bright field of your samples. An additional Osram 100 W Epi fluorescence illumination unit is available to you as a fluorescence microscope (OCM 165) for perfect illumination and stimulation of your fluorescence samples
- **1** A special Abbe N.A. 0.3 condenser with aperture diaphragm and large working distance of 72 mm guarantees the very best working practise in the bright field and with fluorescence applications
- As standard, the OCM range is fitted with a trinocular eyepiece tube
- The mechanical stage including specimen holder (Ø 118 mm) means that you can work quickly and effectively. Further brackets for petri dishes are included with delivery or available as accessories
- Further options such as, for example, a selection of eyepieces, objectives, specimen holders and other phase contrast units can be integrated as accessories
- A protective dust cover as well as multi-lingual user instructions are included in the scope of delivery

Scope of application

- Research and breeding of cell cultures and tissue cultures

Applications/Samples

- Particularly for viewing samples in culture vessels (flasks, petri dishes, microtitre plates), translucent, thin, low-contrast, challenging samples (e.g. immuno fluorescence, FISH, DAPI staining, etc.)

Technical data

- Infinity optical system
- **2** Quintuple nosepiece
- Siedentopf-head, 45° tilt
- Dioptre compensation on both sides

OCM 161


- Overall dimensions W×D×H 304×599×530 mm
- Net weight approx. 13,5 kg

OCM 165

- Overall dimensions W×D×H 304×782×530 mm
- Net weight approx. 21 kg

Please refer to the website for further details as well as information on the full range of model features

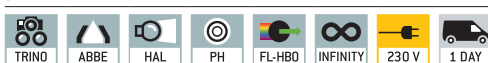STANDARD



Model	Standard configuration					Price excl. of VAT ex works €
	Tube	Eyepiece	Objective quality	Objectives	Illumination	
KERN						
OCM 161	Trinocular	HWF 10×/Ø 22 mm	Infinity Plan	LWD10×/LWD20×/ LWD40×/ LWD20×PH	6 V/30 W Halogen (transmitted)	3180,-
OCM 165	Trinocular	HWF 10×/Ø 22 mm	Infinity Plan		6 V/30 W Halogen + 100 W Epi Fluorescence (B/G)	7790,-




KERN ORF

Refractive index measurement for hospitals, medical training institutions or in sports medicine

! ▲ Also available as digital models or ABBE refractometers, we are happy to advise you

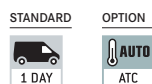
Features

- The KERN ORA refractometers are universal, maintenance-free analogue handheld refractometers
- The main scope of applications is: Hospitals, Doctor's surgeries/Physicians, Medical training institutions, Nursing homes, Sports medicine (doping test)
- Particularly the ORA 2PA is suitable for the measurement of the specific gravity (sg) in urine, the quantity of serum (serumproteine) in urine (doping control among athletes), and the refractive index. The ORA 5PB has been developed specially for the veterinary sector
- The handy and robust design allows the easy, efficient and sustainable use in everyday life.
- Manually calculated conversions and errors of the user are avoided by multiple selectable scales
- These scales are especially developed, exactly calculated and checked. They are also characterized by their thin and clear lines
- The optical system and the prism cover are made of special material which allows a low-tolerance measuring
- All ORA models are equipped with an eyepiece for easy and smooth setting for many different diopter strengths
- The models marked with "ATC" have an automatic temperature compensation which enables accurate measurement at different ambient temperatures
-  The following accessory-parts are included:
 - Storage box
 - Calibration liquid
 - Pipette
 - Small screwdriver
 - Cleaning tissue
- Further accessories are optionally available


Technical data

- Die-cast housing of copper-aluminium alloy, chrome coated
- Measurement temperature without ATC: 20 °C
- Measurement temperature range with ATC: 10 °C/30 °C
- Dimensions of the box W×D×H 205×75×55 mm
- Packing dimensions W×D×H 215×94×65 mm
- Product length: approx. 130–200 mm (depending on the model)
- Net weight approx. 145 g

Please refer to the website for further details as well as information on the full range of model features




Model	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
KERN					
ORA 4RR	Refractive index	1,440 – 1,520 nD	0,001 nD		85,-
ORA 2PA	Serum protein Urine (spec. gravity) Refractive index	0 – 12 g/dl 1,000 – 1,050 sgU 1,3330 – 1,3600 nD	0,2 g/dl 0,002 sgU 0,0005 nD	✓	85,-
ORA 32BA	Brix	0 – 32 %	0,2 %	✓	85,-
ORA 5PB	Serum protein Urine (s. g. dog) Urine (s. g. cat)	2 – 14 g/dl 1,000 – 1,060 sgU 1,000 – 1,060 sgU	0,1 g/dl 0,001 sgU 0,001 sgU		80,-

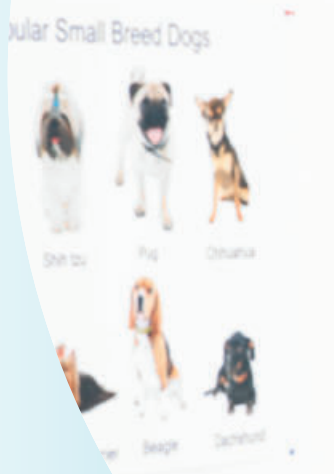

Animal scales for weighing various animals such as rabbits, cats, dogs, snakes, sheep, horses, etc.

Many useful functions assist the vet, animal keeper or pet owner in their daily work with the animal. Thanks to a special animal weighing program, a steady weight will be displayed, even if the animal moves during the weighing process. You can use the practical HOLD function, which means that the weight obtained will continue to be displayed even after the animal has been removed from the scale. This means that you have time to take care of the animal first.

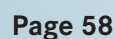KERN – for the love of animals!





Practical assistant for the mobile weighing of heavy pets, animal feed etc. up to 1 t!

Features

- Thanks to the animal weighing function, it can also be used by vets. Compact – so it can be stored in the vet's case and is ideally suited for heavy pets
- Data hold function:
For easy reading of the weight, the display can be "frozen" in different ways. Either automatically when the weighing value remains unchanged or manually by pressing the Hold key
- Varied options for hanging the items, thanks to the stainless steel hook and stainless steel snap links included as standard (not with models HCB with [Max] ≥ 0.5 t)
- Peak load display (peak hold), measuring frequency 5 Hz
- [Max] ≤ 200 kg: Snap link (stainless steel) with safety catch standard
- [Max] ≤ 200 kg: Hook (stainless steel) standard

HCB

- Ready for use: Batteries standard, 3×1,5 V AA, operating time approx. 300 h
- **1** HCB with [Max] ≥ 0,5 t: Reinforced versions for heavy loads. Snap link and hook not included

HCN

- Ready for use: Batteries included 2×1,5 V AA, operating time approx. 200 h
- **2** Housing stainless steel, IP65: Protected against dust and water splashes
- **3** Second display on the rear of the balance

Technical data

- LCD display, digit height 12 mm
- Hole diameter of load support, model with [Max] ≤ 200 kg: approx. 14 mm
[Max] > 200 kg: approx. 16 mm

Accessories for models with [Max] ≥ 100 kg

- Hook with safety catch, cast steel, galvanised and lacquered, revolving. Scope of supplies: 2 shackles, 1 lacquered screw bolt, 1 hook, KERN YHA-06, € 125,-

Accessories for models with [Max] ≤ 200 kg

- Hook (stainless steel), aperture approx. 25 mm, KERN HCB-A02, € 30,-
- Snap link (stainless steel) with safety catch, aperture approx. 15 mm, can be reordered, KERN HCB-A01, € 33,-

Accessories for models with [Max] ≥ 0,5 t

- High-strength shackles, hot-dipped galvanised cast steel bracket, curved edge. Scope of supplies: 2 shackles with lacquered screw bolts, suitable for the following models: KERN HCB with [Max] ≥ 0,5 t, KERN YSC-03, € 50,-

STANDARD



Model	Weighing range	Readability	LCD display, digit height	Dimensions housing W×D×H mm	Height with snap link and hook mm	Net weight approx. kg	Price excl. of VAT ex works €
KERN	[Max] kg	[d] g	mm				
HCB 20K10	20	10	12	80×45×150	290	0,60	105,-
HCB 50K20	50	20	12	80×45×150	290	0,60	115,-
HCB 50K100	50	100	12	80×45×150	290	0,60	90,-
HCB 99K50	99	50	12	80×45×150	290	0,60	100,-
HCB 100K200	100	200	12	80×45×150	290	0,35	99,-
HCB 200K100	200	100	12	80×45×150	290	0,60	115,-
HCB 200K500	200	500	12	80×45×150	290	0,60	95,-
HCB 0.5T-3	500	1000	12	80×42×152	-	0,35	190,-
HCB 1T-3	1000	2000	12	80×45×150	-	0,35	250,-
stainless steel							
HCN 50K100IP	50	100	12	96×48×179	315	0,8	130,-
HCN 100K200IP	100	200	12	96×48×179	315	0,8	130,-
HCN 200K500IP	200	500	12	96×48×179	315	0,8	130,-



All-purpose small pet scale for mobile and stationary use

Features

- Ideal table scale for weighing small animals or for the laboratory
- Particularly flat design
- Plastic weighing tray, removable
- Easy and hygienic cleaning
- Tare function facilitates formulation work
- Simple and convenient 2-key operation

Technical data

- LCD display, digit height 15 mm
- Dimensions of weighing plate:
 - A** Ø 82 mm
 - B** Ø 105 mm
 - C** Ø 150 mm, see larger picture
- Batteries included, 9 V block respectively 2×1,5 V AA
- Socket for mains adapter included

Accessories

- External universal mains adapter, with universal input and optional input socket adapters for EU, GB, USA, not included, KERN YKA-03N, € 32,-

STANDARD



OPTION



Model	Weighing range [Max] g	Readability [d] g	Weighing plate	Dimensions housing W×D×H mm	Net weight approx. kg	Price excl. of VAT ex works €
KERN						
EMB 200-3	200	0,001	A	170×240×54	0,85	210,-
EMB 200-2	200	0,01	B	170×240×39	0,50	125,-
EMB 600-2	600	0,01	B	170×240×39	0,50	160,-
EMB 2000-2	2000	0,01	C	170×240×54	0,90	205,-
EMB 500-1	500	0,1	C	170×240×39	0,55	55,-
EMB 1200-1	1200	0,1	C	170×240×39	0,55	78,-
EMB 6000-1	6000	0,1	C	170×240×39	0,55	160,-
EMB 2200-0	2200	1	C	170×240×39	0,55	57,-
EMB 5.2K1	5200	1	C	170×240×39	0,55	60,-
EMB 5.2K5	5200	5	C	170×240×39	0,55	55,-



KERN FOB-NS



KERN FOB-NL

Compact stainless steel small pet scales with high degree of protection against dust and water splashes

	KERN FOB-NS	KERN FOB-NL
Features	<ul style="list-style-type: none"> Stainless steel design of the housing and weighing plate. Its smooth surfaces make it simple to clean High mobility: thanks to compact, flat construction, it is suitable for use in several locations, secure and non-slip positioning through rubber feet Your support in an HACCP-compliant quality system Protective working cover included with delivery 	
	<ul style="list-style-type: none"> IP65: Protected against dust and water splashes (only when using battery) Particularly flat design 	<ul style="list-style-type: none"> IP67: Protected against dust and water splashes (only when using battery)
Technical data	<ul style="list-style-type: none"> Large LCD display, digit height 20 mm Overall dimensions W×D×H 170×150×40 mm Ready for use: battery 9 V Block standard, operating time up to 24 h Permissible ambient temperature 10 °C/35 °C 	<ul style="list-style-type: none"> Large backlit LCD display, digit height 25 mm Overall dimensions W×D×H 285×255×90 mm Ready for use: battery 4×1,5 V AA standard, operating time up to 66 h without backlight Permissible ambient temperature 5 °C/35 °C
Accessories	<ul style="list-style-type: none"> Protective working cover over keyboard and housing. Scope of delivery 5 pcs., KERN FOB-A12S05, € 40,- Mains adapter external, KERN FOB-A09, € 35,- 	<ul style="list-style-type: none"> Protective working cover over keyboard and housing. Scope of delivery 5 pcs., KERN FOB-A13S05, € 40,- Mains adapter external, not included, can be reordered, KERN FOB-A01, € 35,-

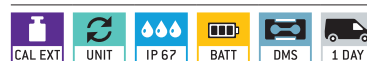
STANDARD



OPTION



STANDARD



OPTION



Model	Weighing range	Readability	Net weight approx.	Weighing plate W×D	Price excl. of VAT ex works €
KERN	[Max] kg	[d] g	kg	mm	
FOB 0.5K-4NS	0,5	0,1	0,80	150×120	105,-
FOB 5K-3NS	5	1	0,80	150×120	105,-
*Dual-range balance switches automatically to the next largest weighing range [Max] and readability [d].					
FOB 3K-4NL	3	0,2	3,80	252×200	265,-
FOB 7K-4NL*	5 7,5	0,5 1	3,80	252×200	230,-
FOB 10K-3NL*	8 15	1 2	3,80	252×200	230,-
FOB 30K-3NL*	16 30	2 5	3,80	252×200	230,-



KERN EOB

KERN EOE



Compact veterinary scale for diverse applications – also with XL platform

Features

- Platform scale to weigh animals
- Stainless steel (EOB) respectively lacquered steel (EOE) weighing plate, which makes cleaning easy and hygienic
- Animal weighing program: When the animal moves during the weighing process, a stable weight is calculated using an average weight and this is then “frozen”. This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- Simple and convenient 4-key operation
- Integrated AUTO-OFF function to preserve the batteries
- Wall mount for the display device
- Protective working cover included with delivery

Technical data

- Large LCD display, digit height 25 mm
- Dimensions of display device W×D×H 235×114×51 mm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 60 h (EOB), 100 h (EOE)
- Mains adapter external, standard



Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A04BS05, **€ 40,-**
- **1** Stand to elevate display device (for models with weighing plate sizes 315×305 mm), height of stand approx. 480 mm, KERN EOB-A01N, **€ 50,-**
- **2** Stand, for all models EOB, EOE, height of stand approx. 1000 mm, KERN EOB-A02B, **€ 80,-**
- KERN EOE: Non-slip rubber mat standard, W×D×H 945×505×5 mm, KERN EOE-A01, **€ 45,-**

STANDARD




Model	Weighing range	Readability	Weighing plate	Cable length approx.	Net weight approx.	Price excl. of VAT ex works €
	[Max]	[d]	mm	cm	kg	
KERN	kg	g				
EOB 15K5	15	5	315×305×57	200	4	152,-
EOB 35K10	35	10	315×305×57	180	4	152,-
EOB 60K20	60	20	315×305×57	220	4	152,-
EOB 60K20L	60	20	550×550×58	270	14	330,-
EOE 150K50XL	150	50	950×500×58	270	18	310,-
EOE 300K100XL	300	100	950×500×58	270	18	310,-

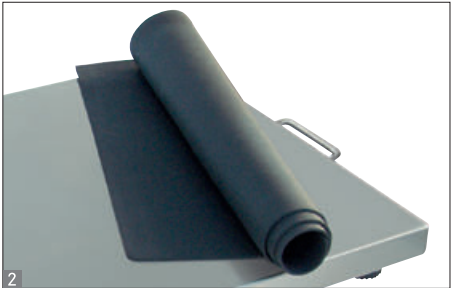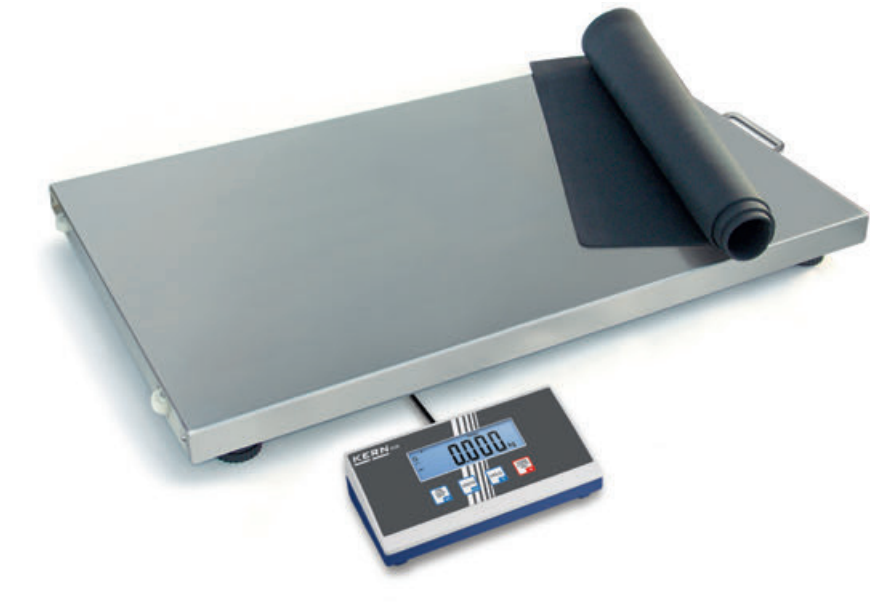

Affordable veterinary model which is suitable for stationary or mobile use


Features

- Robust platform scale to weigh large animals
- Weighing plate made out of stainless steel, for easy and hygienic cleaning
- The low platform height facilitates access
- Animal weighing program: When the animal moves during the weighing process, a stable average value is obtained by pressing the HOLD key. This value remains “frozen” until the HOLD key is pressed again. This leaves sufficient time to look after the small patients and then read the weight in peace
- Simple and convenient 4-key operation
- Integrated AUTO-OFF function to preserve the batteries, adjustable
-  The scale can be easily transported using 2 rollers and a handle and does not require much storage space
- Non-slip rubber mat standard
- Wall mount for display device, standard
- Protective working cover included with delivery

Technical data


- Large LCD display, digit height 25 mm
- Dimensions of display device W×D×H 235×114×51 mm
- Dimensions of weighing plate, stainless steel W×D×H 950×500×58 mm
- Battery operation possible, 4×1,5 V AA, not included
- Mains adapter external, standard


Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A04BS05, **€ 40,-**
- Stand to elevate display device, height of stand approx. 1000 mm, KERN EOB-A02B, **€ 80,-**
-  Non-slip rubber mat, W×D×H 945×505×5 mm, KERN EOE-A01, **€ 45,-**

STANDARD



Model	Weighing range	Readability	Cable length approx.	Net weight approx.	Price excl. of VAT ex works €
KERN	[Max] kg	[d] g	cm	kg	
EOS 150K50XL	150	50	270	17	395,-
EOS 300K100XL	300	100	270	17	395,- 

 Price reduction